


- 4-way stretch antimicrobial fabric surface.
- Continuous layer of 3/8" Akton® polymer.
- Soft viscofoam mid-layer.
- High-density foam base.
- Fluid proof, welded seams.
- 100% no-stitch hook-and-loop attachment.
- Proprietary hidden vent system.
- Easy to clean waterfall construction.
- Available in multiple thicknesses.
- Armboards available.
- Custom pads are available. Call for details.

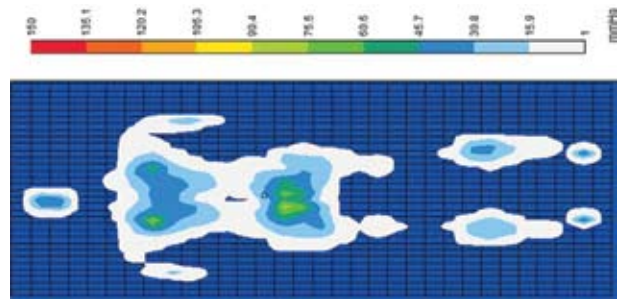

3" Performance Pad  
 5'11" Male Subject, 175 lbs. Average Pressure (mmHg): 18.72 Room Temp: 68 degrees Fahrenheit  
 Equipment: FSA Pressure Map System Range: 1 – 150 mmHg

Introducing our **Performance** and **Performance PLUS Series** of table pads that will help you achieve your zero tolerance goals for HAPUs (Hospital Acquired Pressure Ulcers) and fall well within your budget.


Our pads are clinically engineered to give the maximum in decubiti prevention, infection control, and stability for the ultimate in patient safety. With two styles and multiple thicknesses, we have the perfect mattress for any budgetary requirement.


## Action® Performance Series O.R. Table Pads

<b>55120</b>	2" Fabric Top	Fits Amsco 3080/3085
<b>55130</b>	2" Fabric Top	Fits Steris C-Max
<b>55310</b>	2" Fabric Top	Fits Skytron 3500
<b>55320</b>	2" Fabric Top	Fits Skytron 5001/ 6001/ 6002B/ 6500/ 6600
<b>55330</b>	2" Fabric Top	Fits Skytron 6700/ 6701
<b>55370</b>	2" Fabric Top	Fits Skytron 3501
<b>57120</b>	3" Fabric Top	Fits Amsco 3080/3085
<b>57130</b>	3" Fabric Top	Fits Steris C-Max
<b>57310</b>	3" Fabric Top	Fits Skytron 3500
<b>57320</b>	3" Fabric Top	Fits Skytron 5001/ 6001/ 6002B/ 6500/ 6600
<b>57330</b>	3" Fabric Top	Fits Skytron 6700/ 6701
<b>57370</b>	3" Fabric Top	Fits Skytron 3501
<b>59100</b>	2" Fabric Top	Universal Armboard
<b>59600</b>	3" Fabric Top	Universal Armboard




Durable, welded, midline seam and proprietary hidden vent structure enhance infection control for maximum pad life.

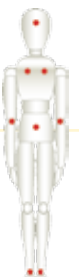

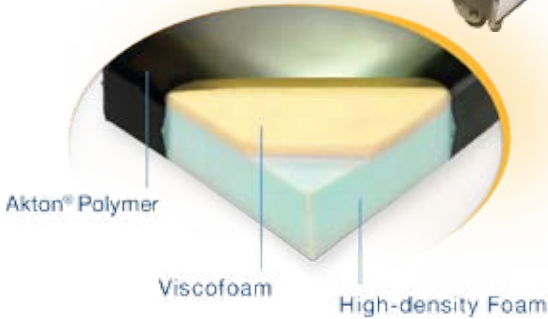

- Shear Smart® surface of exposed 3/8" Akton® polymer.
- Multi-directional stretch film covering on surface.
- Soft viscofoam mid-layer.
- High-density foam base.
- Tissue equivalent.
- Does not support bacterial growth.
- Fluid proof, welded seams.
- 100% no-stitch hook-and-loop attachment.
- Proprietary hidden vent system.
- Easy to clean waterfall construction.
- Available in multiple thicknesses.
- Armboards available.
- Custom pads are available. Call for details.

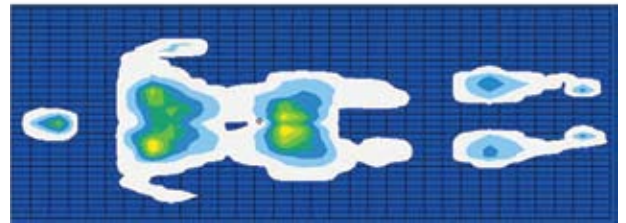

3" Performance Plus Pad  
 5' 11" Male Subject, 175 lbs. Average Pressure (mmHg): 22.33 Room Temp: 69 degrees Fahrenheit  
 Equipment: FSA Pressure Map System Range: 1 – 150 mmHg

## Action® Performance Plus Series O.R. Table Pads

<b>50120</b>	2" Polymer Top	Fits Amsco 3080/3085
<b>50130</b>	2" Polymer Top	Fits Steris C-Max
<b>50310</b>	2" Polymer Top	Fits Skytron 3500
<b>50320</b>	2" Polymer Top	Fits Skytron 5001/ 6001/ 6002B/ 6500/ 6600
<b>50330</b>	2" Polymer Top	Fits Skytron 6700/ 6701
<b>50370</b>	2" Polymer Top	Fits Skytron 3501
<b>52120</b>	3" Polymer Top	Fits Amsco 3080/3085
<b>52130</b>	3" Polymer Top	Fits Steris C-Max
<b>52310</b>	3" Polymer Top	Fits Skytron 3500
<b>52320</b>	3" Polymer Top	Fits Skytron 5001/ 6001/ 6002B/ 6500/ 6600
<b>52330</b>	3" Polymer Top	Fits Skytron 6700/ 6701
<b>52370</b>	3" Polymer Top	Fits Skytron 3501
<b>54100</b>	2" Polymer Top	Universal Armboard
<b>54600</b>	3" Polymer Top	Universal Armboard



AKTON polymer in a proprietary soft film that reduces shear effects for optimal patient safety.



### Resurfacing Kits

Quick and easy resurfacing of minor cuts or abrasions on any polymer surface.

- KIT002** 4 x 4" 10 x 10 cm
- KIT003** 36 x 4" 91 x 10 cm

- VALUE product essential for any operating room.

