


Making a Material Difference[®] with clinically engineered O.R. products
Action Products, Inc.

ACTION[®]

Action[®] Support Surfaces

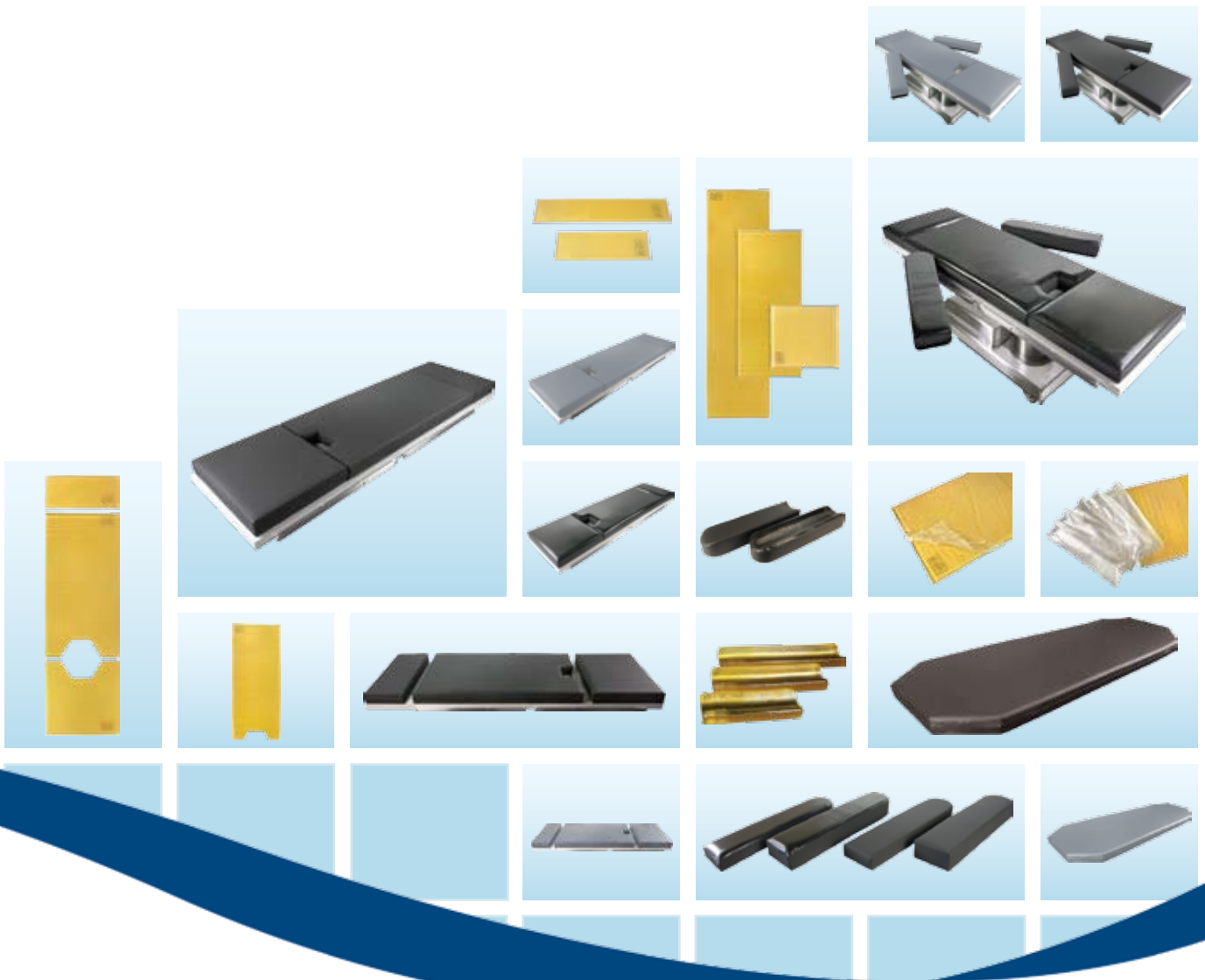

www.actionproducts.com

Making a Material Difference®

The Action® brand is synonymous with clinically engineered quality. Ours is a name you can trust for head-to-toe pressure protection, in and out of the Operating Room. We have designed a variety of support surface options to fit your facility's individual needs. You choose what balance of clinical effectiveness and cost efficiency is right for you. The benefits you can expect from every Action® support surface include:

Infection Control Design

- Bacteria / Microbe resistant
- Midline seams avoid pooling liquids
- Fluid/waterproof, welded seams
- 100% no-stitch construction
- 100% no-stitch hook and loop attachment
- Proprietary hidden vent system minimizes cross contamination



Durable, welded midline seam and proprietary hidden vent structure enhances infection control and extends pad life.



Akton® polymer with a proprietary soft film that reduces shear effects for optimal patient safety.



Why Silver Fabric? Silver has been used for centuries to kill bacteria. A **silver fabric upgrade is available on all Action® support surfaces.** This fabric has active silver ions adhered to a four layer polyurethane coating. This surface is proactive in eliminating bacteria and odors, while reducing bio-burdens within facilities like yours.

- 4-way stretch fabric for shear reduction
- Infection Control Design
- High-density foam base for support
- Imaging compatible
- Meets ergonomic weight limits for hospital staff
- Latex, Silicone and Plasticizer Free
- Static resistive
- CAL133 compliant
- Easy to clean
- Waterfall construction
- Available in multiple thicknesses
- Armboard pads available
- Large selection of models
- Custom pads are available- Call for details



Surface Benefits



Akton® Polymer

- Best shearing qualities
- 360° multi-directional stretch
- Enhances patient stability
- Redistributes heat
- Redistributes pressure
- Repairable
- Bacteria / Microbe resistant



4-way Stretch Fabric

- Use for heavy frames
- Familiar stretch surface
- Bacteria / Microbe resistant solution in fibers of fabric




Silver Fabric

- Proactive infection control
- Odor control
- Bacteria / Microbe resistant silver ion film adhered to top surface
- 4-way stretch fabric


Visual Comparison Chart

4-way Stretch Fabric
 Akton® Polymer
 Viscofoam
 High-density Foam

Performance

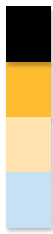

- Premium pressure management
- Great for extended cases
- Frame friendly
- Coordinating flat and contoured armboard pads
- Minimal pressure change over time
- Available in Silver fabric

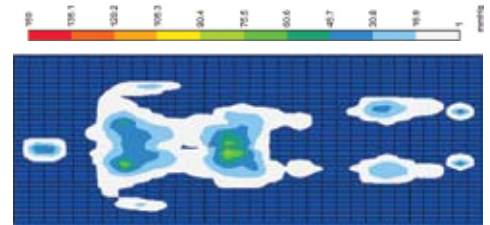


3" Performance Pad • 5'11" Male Subject, 175 lbs.
Average Pressure (mmHg): 18.72 • Room Temp: 68 °F

Performance Plus



- Premium pressure and shear management
- 360° multidirectional stretch surface
- Great for extended cases
- Coordinating flat and contoured armboard pads
- Minimal pressure change over time
- Repairable surface

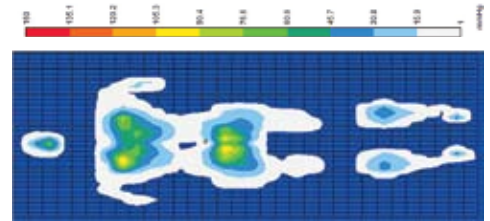


3" Performance Plus Pad • 5'11" Male Subject, 175 lbs.
Average Pressure (mmHg): 22.23 • Room Temp: 69 °F

Original FT



- Action® Value Line Product
- Good pressure management
- Good for extended cases
- Frame friendly
- Coordinating flat armboard pads
- Available in Silver fabric

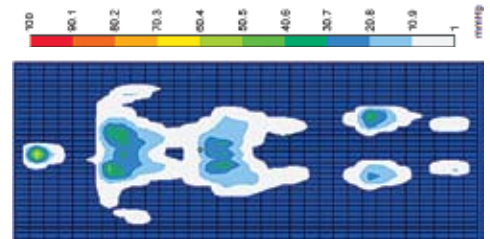


3" Original FT Pad • 5'7" Male Subject, 170 lbs.
Average Pressure (mmHg): 20.58 • Room Temp: 70 °F

Original



- Action® Value Line Product
- 360° multi-directional stretch surface
- Premium shear management
- Good pressure management
- Good for extended cases
- Coordinating flat armboard pads
- Repairable surface

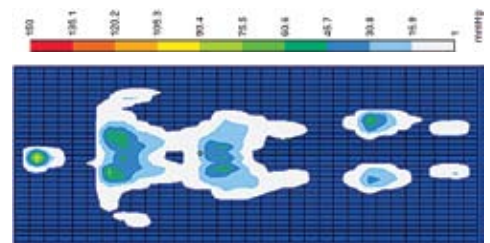


3" Original Pad • 5'11" Male Subject, 180 lbs.
Average Pressure (mmHg): 19.4 • Room Temp: 70 °F

Angel





- Premium infection control pad
- Good pressure management
- Good for short-term cases
- Coordinating flat and contoured armboard pads
- Available in Stretcher Pads
- Available in Silver fabric

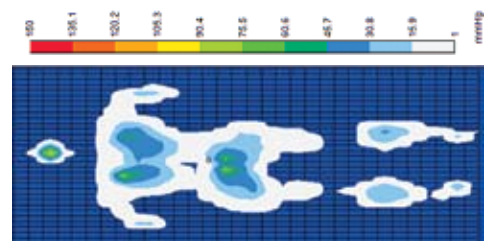

3" Angel Pad • 5'11" Male Subject, 180 lbs.
Average Pressure (mmHg): 16.34 • Room Temp: 70 °F

Action® Performance and Performance Plus Table Pads

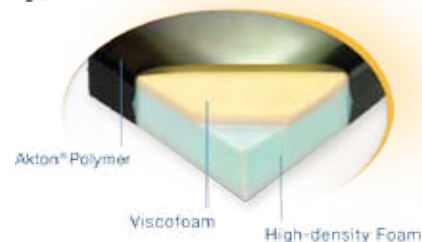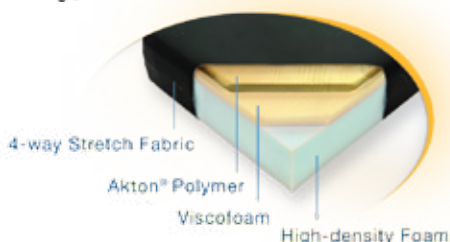

Action® Performance Series

Head-to-toe pressure management, perfect for extended cases

- 4-way stretch fabric top provides pressure and shear management
- Continuous layer of exclusive Akton® viscoelastic polymer maximizes pressure redistribution and prevents capillary closure
- Soft viscofoam mid-layer enhances patient envelopment and comfort
- Superior infection control design
- High-density foam base maximizes stability

Action® Performance Plus Series

The Performance Plus Series provides all of the pressure management features of the Performance Series- **PLUS** the added benefit of:

- Exposed Akton® polymer surface provides a shear benefit that is unparalleled in the industry
- Exclusive soft, 360° multi-directional stretch film cover preserves skin integrity
- Redistribution of heat to prevent “hot spots” which can become a catalyst for decubiti
- Repairable polymer top
- Clinical studies recommend a viscoelastic polymer overlay atop foam pads for effective pressure redistribution*



Models	Performance	Performance	Performance Silver	Performance Silver	Performance Plus	Performance Plus
	2 3/8" 3-pc. Set	3 3/8" 3-pc. Set	2 3/8" 3-pc. Set	3 3/8" 3-pc. Set	2 3/8" 3-pc. Set	3 3/8" 3-pc. Set
Steris Amsco 3080/3085	55120	57120	55120DS	57120DS	50120	52120
Steris C-Max	55130	57130	55130DS	57130DS	50130	52130
Skytron 3500	55310	57310	55310DS	57310DS	50310	52310
Skytron 5001, 6001, 6002, 6500 & 6600	55320	57320	55320DS	57320DS	50320	52320
Skytron 6700, 6701	55330	57330	55330DS	57330DS	50330	52330
Skytron 3501, 3600, 3600B	55370	57370	553710DS	57370DS	50370	52370

Don't see what you are looking for.... Call customer service for additional sizes or email service@actionproducts.com.

Action® Original FT and Original Series Table Pads



Value Line - Essential O.R. Products

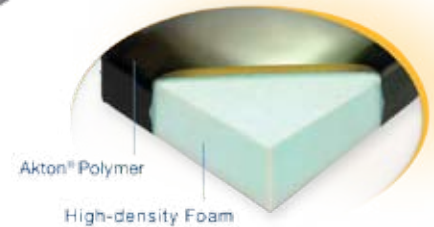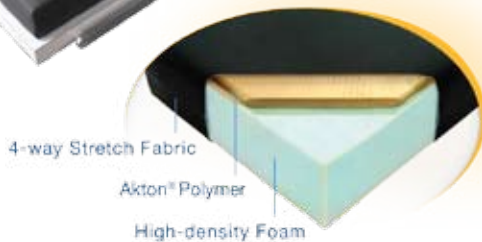

Action® Original FT Series

Good pressure management with our frame friendly, infection control design

- Value pad designed for a budget-driven facility
- 4- way stretch fabric top provides pressure and shear management
- Continuous layer of exclusive Akton® viscoelastic polymer maximizes pressure redistribution and prevents capillary closure
- Superior infection control design
- High-density foam base maximizes stability

Action® Original Series

Flagship pad with upgraded infection control design

- Exposed Akton® polymer surface providing shear benefit that is unparalleled in the industry
- Exclusive soft, 360° multi-directional stretch film cover preserves skin integrity
- Redistribution of heat to prevent “hot spots” which can become a catalyst for decubiti
- Repairable polymer top
- Superior infection control design
- High-density foam base maximizes stability
- Clinical studies recommend a viscoelastic polymer overlay atop foam pad for effective pressure redistribution*

Models	2 3/8" Original FT-	3 3/8" Original FT-	2 3/8" Original FT Silver	3 3/8" Original FT Silver	2 3/8" Original	3 3/8" Original
	3-pc. Set	3-pc. Set	3-pc. Set	3-pc. Set	3-pc. Set	3-pc. Set
Steris Amsco 3080/3085	75120	77120	75120DS	77120DS	70120	72120
Steris C-Max	75130	77130	75130DS	77130DS	70130	72130
Skytron 3500	75310	77310	75310DS	77310DS	70310	72310
Skytron 5001, 6001, 6002, 6500 & 6600	75320	77320	75320DS	77320DS	70320	72320
Skytron 6700, 6701	75330	77330	75330DS	77330DS	70330	72330
Skytron 3501, 3600, 3600B	75370	77370	75370DS	77370DS	70370	72370


Don't see what you are looking for.... Call customer service for additional sizes, or email service@actionproducts.com.

Action® Angel O.R. Table Pads / Overlays

Action® Angel Operating Room Table Pads

Outstanding replacement pad for today's fiscally sensitive healthcare environment. Action® Angel series features superior infection control construction and good pressure and shear management performance. This is a Body Smart® and budget friendly alternative to the "basic table pad". It was engineered specifically for the dynamic, high turnover demands of the busy surgical suite without sacrificing patient care. Perfect for:

- ASCs
- Higher turnover suites and procedures
- Use with existing Action® Overlays
- Comfort for non-anesthetized patient procedures
- Affordable upgrade alternative to the "basic" pad

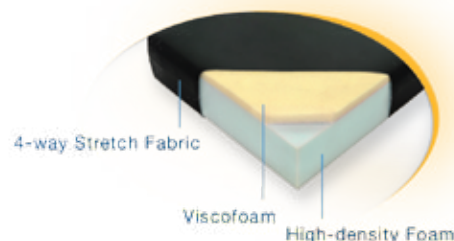

Models	2" 3-pc. Set	3" 3-pc. Set	4" 3-pc. Set	2" Silver 3-pc. Set	3" Silver 3-pc. Set	4" Silver 3-pc. Set
Steris Amsco 3080/3085	60120	62120	65120	60120DS	62120DS	65120DS
Steris C-Max	60130	62130	65130	60130DS	62130DS	65130DS
Skytron 3500	60310	62310	65310	60310DS	62310DS	65310DS
Skytron 5001, 6001, 6002, 6500 & 6600	60320	62320	65320	60320DS	62320DS	65320DS
Skytron 6700, 6701	60330	62330	65330	60330DS	62330DS	65330DS
Skytron 3501, 3600, 3600B	60370	62370	65370	60370DS	62370DS	65370DS

Don't see what you are looking for.... Call customer service for additional sizes, or email service@actionproducts.com.

Action® O.R. Overlays

The Action® O.R. Overlay is a versatile pad used for a wide range of procedures. It is especially important for use in procedures lasting more than one hour or for patients with skin integrity or circulatory problems. The Akton® polymer pads, with the Shear Smart® 360° multi-directional stretch surface, are suggested for all surgical procedures, including robotics and interventional radiology uses. The Action® Overlay provides pressure, weight and temperature redistribution in a simple, easy to use and clean design. A full length pad or three piece set can extend the life of your standard table pad.



Models	Size	SKU
Action® O.R. Overlay (Small)	20" x 20" x 1/2" (50.8 x 50.8 x 1.3 cm)	40105
Action® O.R. Overlay (Medium)	20" x 46" x 1/2" (117 x 50 x 1.3 cm)	40100
Action® O.R. Overlay (Full Length)	20" x 72" x 1/2" (50 x 183 x 1.3 cm)	40101
Action® O.R. Overlay, Perineal Cutout	20" x 46" x 1/2" with cutout (50.8 x 117 x 1.3 cm)	40103
Action® O.R. Overlay 3-piece Set	See A, B and C	40104
Action® O.R. Overlay Head	20" x 9" x 1/2" (50.8 x 22.9 x 1.3 cm)	40104A
Action® O.R. Overlay Torso	20" x 39 1/2" x 1/2" with perineal cutout (50.8 x 100.3 x 1.3 cm)	40104B
Action® O.R. Overlay Foot	20" x 20" x 1/2" with perineal cutout (50.8 x 50.8 x 1.3 cm)	40104C

Action® Armboard Pads / Overlay Covers

Action® Armboard Pads

All Action® armboard pads are constructed the same as our table and stretcher pads with 4-way stretch fabric, midline seams, no-stitch construction and proprietary hidden vent systems. They are available in a variety of shapes and sizes, and include the hook and loop attachment, to meet your facility's pressure and infection control management needs.

Universal Flat Armboard Pads

Models	Angel			Original FT		Original		Performance			Performance Plus		
	2"	3"	4"	2 3/8"	3 3/8"	2 3/8"	3 3/8"	2 3/8"	3 3/8"	4 3/8"	2 3/8"	3 3/8"	4 3/8"
Flat	64100	64600	67100	79100	79600	74100	74600	59100	59600	59700	54100	54600	54700

Silver fabric has a \$35 upcharge available upon request.


Contoured Armboard Pads

Models	Angel			Performance		Performance Plus	
	2"	3"	4"	2 3/8"	3 3/8"	2 3/8"	3 3/8"
Contoured	64300	64800	67300	59300	59800	54300	54800

Silver fabric has a \$35 upcharge available upon request.

Polymer Armboard Overlay Pads

Models	Size	SKU
Short Armboard Overlay Pad	4" x 15" x 1/4" (10 x 38 x .6 cm)	40300
Long Armboard Overlay Pad	4 1/2" x 24" x 3/8" (11.5 x 61 x .95 cm)	40301
Wide Short Armboard Overlay Pad	15" x 6 1/2" x 1/4" (38 x 16.5 x .6 cm)	40309
Wide Long Armboard Overlay Pad	24" x 6 1/2" x 3/8" (61 x 16.5 x .95 cm)	40308
Small Contoured Armboard Pad	16" x 5 1/2" x 2" (41 x 14 x 5 cm)	40307
Medium Contoured Armboard Pad	20" x 5 3/4" x 2 1/4" (51 x 15 x 6 cm)	40310
Large Contoured Armboard Pad	24" x 5 3/4" x 2 1/4" (61 x 15 x 6 cm)	40306

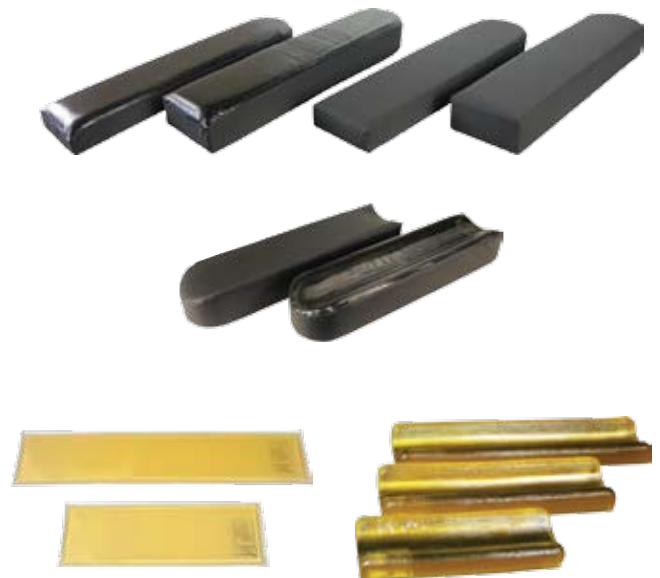

Disposable Overlay Covers Available!

Action® Disposable Overlay Covers aid in post-procedural cleaning and room turnover times and enhances the service life of Action® Overlays. They are made with the same film as our polymer products. With a soft, 360° multi-directional stretch surface, it has great shear benefits as well.



Models	Size	SKU
Action® Disposable Overlay Cover - Pillowcase Design Style	For Models 40100 and 40101- Case of 25	40102
Action® Disposable Overlay Cover - Fitted Sheet Style	For Model 40100 Table Pad- Case of 25	40800
Action® Disposable Overlay Cover - Fitted Sheet Style	For Model 40101 Table Pad- Case of 25	40801
Action® Disposable Overlay Cover - Fitted Sheet Style	For Model 40103 Table Pad- Case of 25	40803
Action® Disposable Overlay Cover - Fitted Sheet Style for 3 Pc	For Model 40104 Table Pad	40804
Action® Disposable Overlay Cover - Fitted Sheet Style	For Model 40105 Table Pad	40805



Action® Angel Stretcher Pads

All patient support surfaces should be safe for the patient and the hospital staff! Action realizes that patients must often spend hours waiting on uncomfortable stretchers that provide zero to minimal pressure or shear relief. Unfortunately, the long wait times can drive the onset of hospital acquired pressure ulcers. Fortunately, Action has developed an affordable solution. Our Angel Stretcher pads offer all the pressure management and infection control benefits of our Angel OR table pads; welded, fluid-proof seams, 4-way stretch, bacteria resistant fabric and a dual-density foam core in either 3" or 4" heights. Consider upgrading and standardizing your support surfaces throughout your facility today!

Models	3" Angel	4" Angel	3" Angel Silver	4" Angel Silver
Universal 24 x 76	67124	68124	67124DS	68124DS
Universal 26 x 76	67126	68126	67126DS	68126DS
Universal 28 x 76	67128	68128	67128DS	68128DS
Universal 30 x 76	67130	68130	67130DS	68130DS

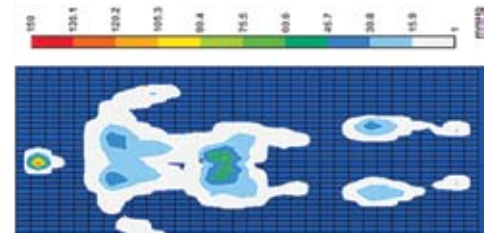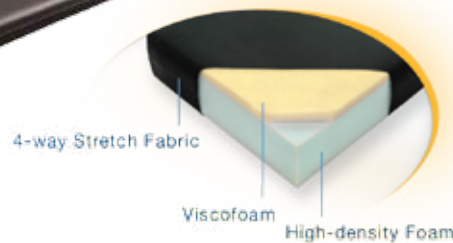

3" Angel Stretcher Pad • 5'11" Male Subject, 180 lbs.
Average Pressure (mmHg): 15.83 • Room Temp: 70 °F

Action Products, Inc.

954 Sweeney Drive
Hagerstown, MD 21740, USA
800.228.7763
301.797.1414
301.733.2073 (Fax)
service@actionproducts.com
www.actionproducts.com
blog.actionproducts.com



Warranty

Action® pads are warranted to be free of defective materials and workmanship at the time when first purchased and are serviceable for a period of two years provided the product is not subject to misuse, negligence, accident or abuse. Specific warranty questions should be directed to Customer Service at 301.797.1414 (Monday - Friday 8:30AM - 5:00PM EST) or your authorized dealer/distributor.

CUSTOM PADS AND POSITIONERS ARE OUR SPECIALTY.

Please call or fax a diagram with your specific requirements. Custom orders cannot be returned. All pads and positioners are imaging compatible and easily disinfected.

TERMS:

Net 30 days, 2% discount on merchandise only for payment within 10 days of invoice date. FOB Hagerstown, MD. Service Charge of 1.5% per month on unpaid balance. Prices may change without notice. Handling charge of \$5.00 per order. \$3.00 drop ship fee may apply.

ACTION, AKTON, SHEAR SMART, BODY SMART and Making A Material Difference are registered trademarks of Action Products, Inc. USA

* AORN Recommended Practices for Positioning the Patient in the Perioperative Setting, 2011, p.337—source from Research in Nursing Health, 1994, 17, p.337, Intraoperative Pressure Sore Prevention: An Analysis of Bedding Materials Thomas Jefferson Study

