



Support Surfaces

O.R. Table Pads & Overlays



Action Products, Inc.
www.actionproducts.com



The Action® brand is a name you can trust for head-to-toe pressure protection, in and out of the Operating Room. We have designed a variety of support surface options to fit your facility's individual needs. You choose what balance of clinical effectiveness and cost efficiency is right for you.



Benefits you can expect from every Action® OR Table Pad:

- 4-way stretch fabric for shear reduction and immersion benefits
- Biocompatible
- PVC free
- Imaging compatible
- Silicone and Plasticizer Free
- Not made with natural rubber latex
- Meets ergonomic weight limits for hospital staff
- Easy to clean
- Waterfall construction
- Available in multiple thicknesses
- Armboard pads available
- Large selection of models

Surface Benefits

4-way Stretch Fabric

- Premium PU-coated waterproof fabric
- Enhances patient immersion
- Fungistatic
- Bacteria / Microbe resistant solution in fibers of fabric
- Use for heavy frames

Akton® Polymer

- Premium shear management
- 360° multi-directional stretch
- Enhances patient stability
- Evenly distributes heat
- Evenly distributes pressure
- Does not support bacterial growth
- Hypoallergenic


Easy To Clean

- Waterproof
- Bacteria / Microbe resistant
- Midline seams avoid pooling liquids
- Welded seams prevent ingress of fluids
- 100% no-stitch construction
- 100% no-stitch hook and loop attachment
- Proprietary hidden vent system helps to minimize cross contamination

Durable, welded midline seam and proprietary hidden vent structure helps to minimize cross contamination and extends pad life.



Welded Fabric Seam and Vent

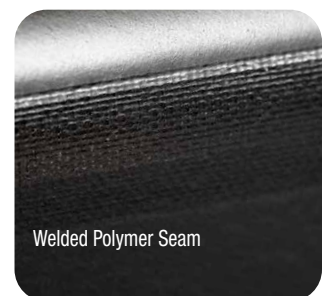

Welded Polymer Seam



Akton® polymer tops have a soft multi-directional film that helps reduce shear effects for patient protection.

AORN 2023 Pressure Injury Prevention Guideline 3.3.1:

“ The researchers found that both the foam and gel mattress and the viscoelastic gel overlay were significantly more effective than the foam mattress in preventing skin changes and pressure injury, but the **viscoelastic gel overlay was the most effective surface for preventing pressure injury.**



Our Performance Series and Original Series table pads have a 3/8" of AKTON polymer built in, so you get a high density foam and overlay in a single pad.

	Best				Better				Good		
	Performance Plus 3-pc. Set		Performance 3-pc. Set		Original 3-pc. Set		Original FT 3-pc. Set		Angel 3-pc. Set		
Pad Surface Type	AKTON Polymer		4-way Stretch Fabric		AKTON Polymer		4-way Stretch Fabric		4-way Stretch Fabric		
Models	2 3/8"	3 3/8"	2 3/8"	3 3/8"	2 3/8"	3 3/8"	2 3/8"	3 3/8"	2"	3"	4"
Berchtold B710, B810	50920	52920	55920	57920	70920	72920	75920	77920	60920	62920	65920
Skytron 3100	50340	52340	55340	57340	70340	72340	75340	77340	60340	62340	65340
Skytron 3500	50310	52310	55310	57310	70310	72310	75310	77310	60310	62310	65310
Skytron 3501, 3600, 3600B, 3502	50370	52370	55370	57370	70370	72370	75370	77370	60370	62370	65370
Skytron 6300, 6302	50350	52350	55350	57350	70350	72350	75350	77350	60350	62350	65350
Skytron 5001, 6001, 6002, 6500 & 6600	50320	52320	55320	57320	70320	72320	75320	77320	60320	62320	65320
Skytron 6700, 6701, 6702	50330	52330	55330	57330	70330	72330	75330	77330	60330	62330	65330
Steris Amsco 1080, 2080	50110	52110	55110	57110	70110	72110	75110	77110	60110	62110	65110
Steris Amsco 3080, 3085	50120	52120	55120	57120	70120	72120	75120	77120	60120	62120	65120
Steris C-Max, 4085, 5085	50140	52140	55140	57140	70140	72140	75140	77140	60140	62140	65140
Berchtold D750/ D850											
Getinge Maquet Meera											
Getinge Maquet Alphamaxx											
Getinge Maquet AlphaStar											
Trumpf 7000dV											

-- Call for item numbers and pricing --
800.228.7763 or 301.797.1414



Have questions about replacement pads for your operating room table?

Email customer service service@actionproducts.com or call 800.228.7763 or 301.797.1414

► Matching Armboards Available - see page 10

The Akton® polymer surface has a soft multi-directional film that helps reduce shear effects for patient protection.



Performance Plus Table Pad shown with other Action® Pads & Positioners

Best

Our Performance Plus Series table pads have three pressure reducing layers with an Akton® polymer surface.

- 1 3/8" layer of **Akton® polymer surface**
- 2 1/2" layer of **viscofoam**
- 3 **high-density foam** base

Performance Plus Series




Top-of-the-line series for pressure and shear management

Premium pressure and shear management benefits of our head-to-toe AKTON polymer surface. Great choice for extended cases.

AKTON Polymer with 360° multi-directional stretch film surface

Viscofoam

High-density Foam

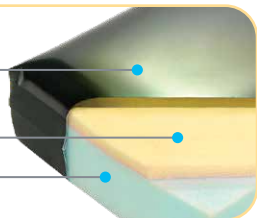

Pressure Redistribution


- 3/8" head-to-toe layer of AKTON viscoelastic polymer surface provides immersion benefits by distributing weight
- Soft viscofoam mid-layer enhances patient envelopment and comfort
- High density foam base provides stability

Temperature Management

- AKTON polymer distributes heat to reduce "hot spots" that are catalysts for pressure injuries and can help reduce heat-related moisture build-up

Shear Management

- Exposed continuous layer of our proprietary AKTON polymer surface provides excellent shear benefits and will not support bacterial growth
- The soft, tissue equivalent polymer surface moves easily with the patient, reducing friction and shear forces

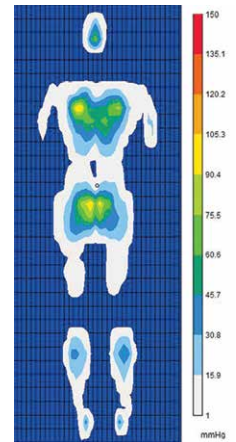


Pressure redistribution of a 3" Performance Plus Pad

► Find the Pad that works with your OR Table Model - see page 3

► Matching Armboards Available - see page 10

Fabric surface
Features a 4-way stretch fabric for increased durability.

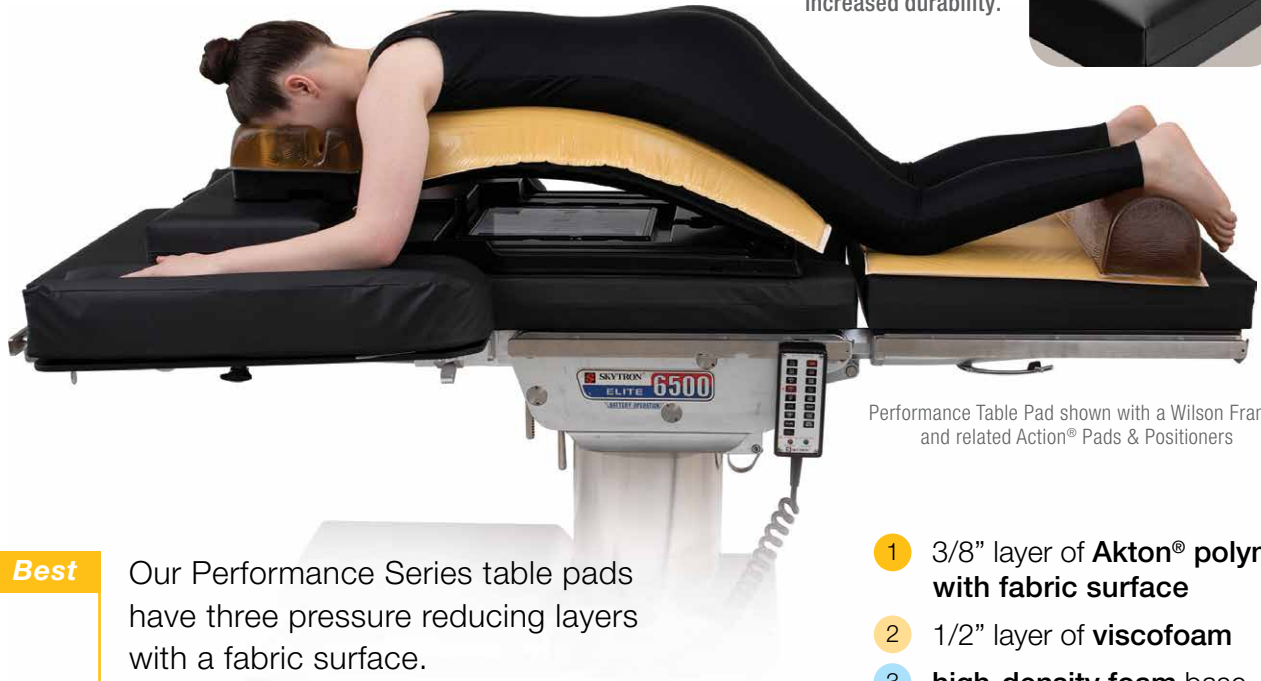

Performance Table Pad shown with a Wilson Frame, and related Action® Pads & Positioners

Best

Our Performance Series table pads have three pressure reducing layers with a fabric surface.

- 1 3/8" layer of **Akton® polymer with fabric surface**
- 2 1/2" layer of **viscofoam**
- 3 **high-density foam base**


Performance Series



Head-to-toe pressure management, great for extended cases and high turnover OR's

Durable 4-way stretch fabric features bacteria/microbe resistant solution in fibers of the fabric.

4-way Stretch Fabric
AKTON Polymer
Viscofoam
High-density Foam

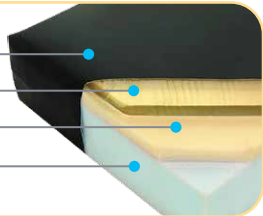

Pressure Redistribution


- 4-way stretch fabric cover compliments pad's ability to immerse the patient; reduces hammocking effect
- Encased 3/8" head-to-toe layer of AKTON viscoelastic polymer provides immersion benefits by distributing weight
- Soft viscofoam mid-layer enhances patient envelopment and comfort
- High density foam base provides stability

Moisture Management

- Vapour release fabric helps to maintain healthy microclimate

Shear Management

- 4-way stretch fabric provides shear management benefits

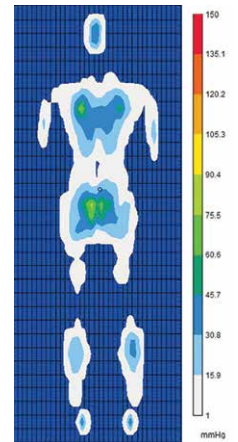

Pressure redistribution of a 3" Performance Pad

► Find the Pad that works with your OR Table Model - see page 3

► Matching Armboards Available - see page 10

The Akton® polymer surface has a soft multi-directional film that helps reduce shear effects for patient protection.



Better

Our Original Series has two pressure reducing layers with an Akton® polymer surface.


- 1 3/8" layer of **Akton® polymer surface**
- 2 **high-density foam base**

Original Series

Premium pressure and shear management benefits of our head-to-toe AKTON polymer surface.

AKTON Polymer with 360° multi-directional stretch film surface

High-density Foam

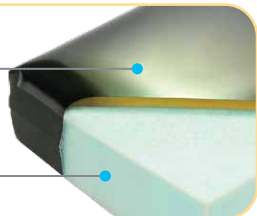

Pressure Redistribution


- 3/8" head-to-toe layer of AKTON viscoelastic polymer surface provides immersion benefits by distributing weight
- High density foam base provides stability

Temperature Management

- AKTON polymer distributes heat to reduce "hot spots" that are catalysts for pressure injuries and can help reduce heat-related moisture build-up

Shear Management

- Exposed continuous layer of our proprietary AKTON polymer surface provides excellent shear benefits and will not support bacterial growth
- The soft, tissue equivalent polymer surface moves easily with the patient; reducing friction and shear forces

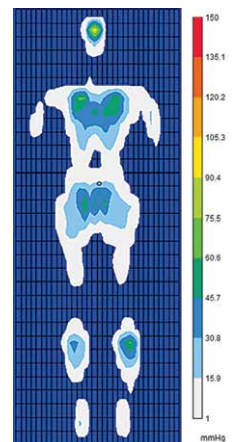

Pressure redistribution of a 3" Original Pad

► Find the Pad that works with your OR Table Model - see page 3

► Matching Armboards Available - see page 10

Fabric Surface
Features a 4-way stretch fabric for increased durability.



Original FT Table Pad shown with a Wilson Frame, and related Action® Pads & Positioners

Better


Our Original FT Series has two pressure reducing layers with a fabric surface.

- 1 3/8" layer of **Akton®** polymer with fabric surface
- 2 high-density foam base

Original FT Series

Durable, budget-friendly pressure management

Durable 4-way stretch fabric features bacteria/microbe resistant solution in fibers of the fabric.

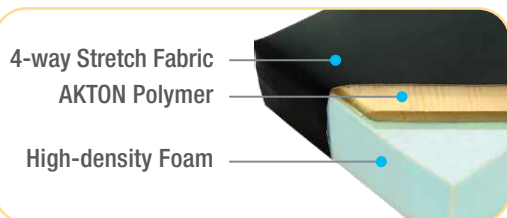

Pressure Redistribution


- 4-way stretch fabric cover compliments pad's ability to immerse the patient; reduces hammocking effect
- Encased 3/8" head-to-toe layer of AKTON viscoelastic polymer provides immersion benefits by distributing weight
- High density foam base provides stability

Moisture Management

- Vapour release fabric helps to maintain healthy microclimate

Shear Management

- 4-way stretch fabric provides shear management benefits

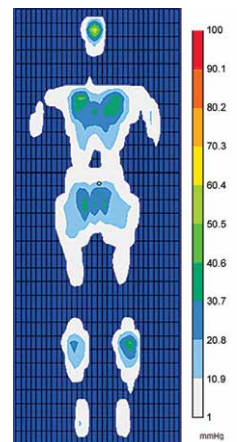

Pressure redistribution of a 3" Original FT Pad

► Find the Pad that works with your OR Table Model - see page 3

► Matching Armboards Available - see page 10



Angel Table Pad shown with our Full Length Overlay and other Action® Pads & Positioners

Good

This is a great budget friendly upgrade of the “basic foam” table pad. It was engineered specifically for the high turnover demands of the busy surgical suite.

- 1 viscofoam with fabric surface
- 2 high-density foam base

Adding an Action® Overlay for longer procedures makes this a very versatile pad

- see page 9

Angel Series

Our budget friendly upgrade of the “basic foam” pad

Replacement Pad for today’s fiscally sensitive healthcare facilities.


Perfect for:

- ASCs
- Higher turnover suites and procedures
- Comfort for non-anesthetized patient procedures

4-way Stretch Fabric

Viscofoam

High-density Foam

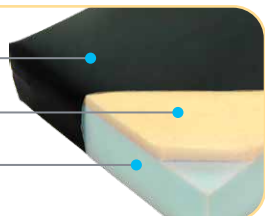

Pressure Redistribution

- 4-way stretch fabric cover compliments pad’s ability to immerse the patient; reduces hammocking effect
- Soft viscofoam layer enhances patient envelopment and comfort
- High density foam base provides stability


Moisture Management

- Vapour release fabric helps to maintain healthy microclimate

Shear Management

- 4-way stretch fabric provides shear management benefits

Combine with an Action® overlay for an improved shear management option for a wider range of patients; like those with fragile skin.

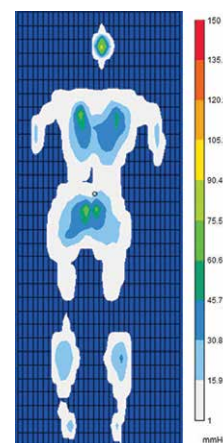

Pressure redistribution of a 3" Angel Pad


► Find the Pad that works with your OR Table Model - see page 3

► Matching Armboards Available - see page 10

Overlays

The Action® O.R. Overlay is a versatile polymer pad used for a wide range of procedures. Its shear reducing qualities combined with the ability to distribute heat and pressure make it a valuable pressure management tool. Used in longer procedures or for patients with fragile skin integrity, the overlay is an affordable upgrade for basic foam table pads that provides enhanced patient protection. The overlays are portable and can be exchanged between rooms as needed and are available in all standard OR table footprints in 1/2" thickness. Disposable covers are available for select overlay models.

40105	Small	20 x 20 x 1/2"	50.8 x 50.8 x 1.3 cm
40100	Medium	46 x 20 x 1/2"	117 x 50 x 1.3 cm
40109	Medium Large	63 x 20 x 1/2"	160 x 50 x 1.3 cm
40101	Full Length	72 x 20 x 1/2"	183 x 50 x 1.3 cm
40103	Perineal Cutout	46 x 20 x 1/2" w/ cutout	117 x 50.8 x 1.3 cm w/ cutout
40104	3-piece Set	See A, B and C	
40104A	Head	9 x 20 x 1/2"	22.9 x 50.8 x 1.3 cm
40104B	Torso w/ perineal cutout	39 1/2 x 20 x 1/2"	100.3 x 50.8 x 1.3 cm
40104C	Foot w/ perineal cutout	20 x 20 x 1/2"	50.8 x 50.8 x 1.3 cm

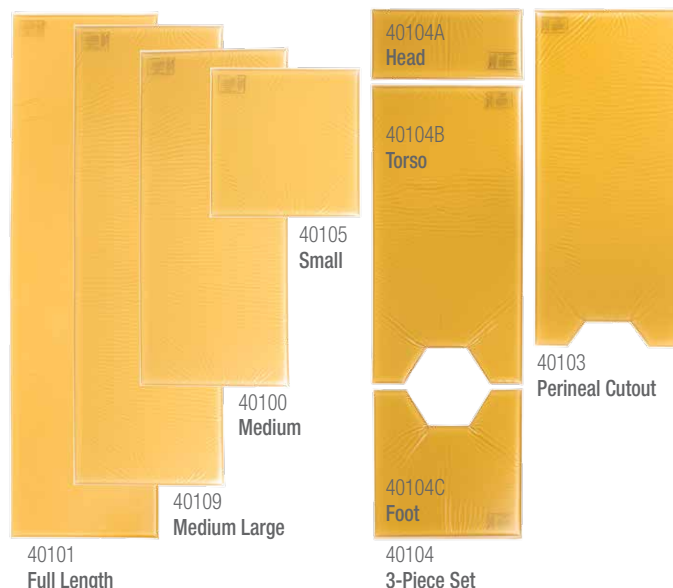


AORN 2023 Guidelines, Pressure Injury Prevention
Selecting and Using Support Surfaces, section 3.2

“Position patients on surfaces that reduce the potential for PI by redistributing pressure and reducing the risk of unintended patient movement... Static surfaces such as viscoelastic polymer (i.e., gel) overlays are designed to reduce shearing and to support the patient's weight without becoming fully compressed under the weight of the patient's body.^{5,46,74} Gel overlays are radiolucent, latex-free, and reusable.⁷⁴



Disposable Overlay Covers

Action's waterproof, non-sterile Disposable Overlay Covers are made of a soft, flexible film and are available in two styles. Designed to help reduce exposure to fluids and to help extend the life of the overlay. The "fitted sheet" style features pockets at both ends of the pad and has a 4" film drape along the length of the overlay. The "pillowcase" style cover fully encases the overlay. This disposable cover will fit the 40100 and 40101 Overlay. *Purchase in cases of 25.*

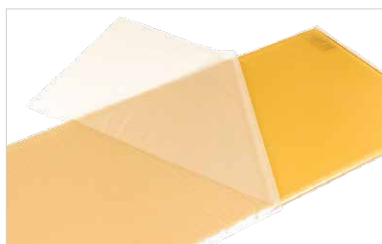

Pillowcase style--full length overlay cover



Fitted Sheet style--torso with cut out

Fitted Sheet Style

40800	Made for overlay model 40100- torso overlay	Case of 25
40801	Made for overlay model 40101- full length pad	Case of 25
40803	Made for overlay model 40103- torso with cut out	Case of 25
40805	Made for small hip pad	Case of 25

Pillowcase Style

40102	Made for table pad models 40100 & 40101	Case of 25
--------------	---	------------


Coordinating Armboard Pads

All Action® armboard pads are constructed the same as our table and stretcher pads with 4-way stretch fabric, midline seams, no-stitch construction and proprietary hidden vent systems. They are available in a variety of shapes and sizes, and include the hook and loop attachment, to meet your facility's pressure management needs.



Models	Flat				Contoured			Incline	
	2 3/8"	3 3/8"	4 3/8"	5 3/8"	2 3/8"	3 3/8"	4 3/8"	2"	3"
Performance Plus	54100	54600	54700	54000	54300	54800	-	N/A	N/A
Performance	59100	59600	59700	59000	59300	59800	59900	59150	59650
Original	74100	74600	N/A	N/A	N/A	N/A	-	N/A	N/A
Original FT	79100	79600	N/A	N/A	N/A	N/A	-	N/A	N/A

Models	Flat			Contoured		
	2"	3"	4"	2"	3"	4"
Angel	64100	64600	67100	64300	64800	67300

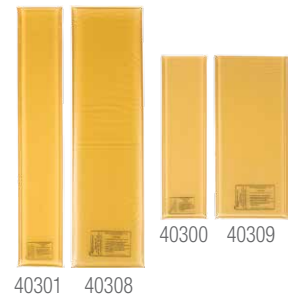

Contoured Armboard Pads

These pads offer excellent full-arm coverage. Contoured design cradles patient's arm, protects bicep, elbows and forearms.

40307	Small	16 x 5 3/4 x 2 1/4"	41 x 14.6 x 5.7 cm
40310	Medium	20 x 5 3/4 x 2 1/4"	51 x 14.6 x 5.7 cm
40306	Large	24 x 5 3/4 x 2 1/4"	61 x 14.6 x 5.7 cm

Armboard Pads

For use on standard armboard table extensions to protect the patient's bicep, elbow and forearm areas.

40300	Short	15 x 4 x 1/4"	38 x 10.2 x .6 cm
40301	Long	24 x 4 x 3/8"	61 x 10.2 x 1 cm
40309	Wide Short	15 x 6 1/2 x 1/4"	38 x 16.5 x .6 cm
40308	Wide Long	24 x 6 1/2 x 3/8"	61 x 16.5 x 1 cm

Stabilizer Straps

The Stabilizer Strap secures and protects the arm.

40315	Short	25 x 2 x 1/4"	63.5 x 5 x .635 cm
40317	Long	36 x 2 x 1/4"	91.4 x 5 x .635 cm



Donut Head Pads

Designed to protect and cradle the patient's head, these pads afford easy access for surgical area of the head and neck. It can be used for a variety of procedures including ENT, plastic, general and ophthalmic.



40201	Adult	8 OD x 3 ID x 1 3/4"	20.5 OD x 7.5 ID x 4.5 cm
40201D	Adult w/ Center Dish	8 OD x 3 ID x 1 3/4"	20.5 OD x 7.5 ID x 4.5 cm
40218	Adult Wide	8 OD x 4 1/2 ID x 1 3/4"	20.5 OD x 11.4 ID x 4.5 cm



Ulnar Forearm Protector

Designed to be wrapped around the arm and secured with the Hook & Loop attachment as the patient's arm is tucked to their side.

40303	18 x 13 x 1/4"	46 x 33 x .6 cm
--------------	----------------	-----------------



Knee Crutch Pads

The Knee Crutch Pad is available in a variety of sizes to protect patient skin integrity.

40402	Small	11 x 9 x 3/8"	27.9 x 22.9 x 1 cm
40403	Medium	17 x 11 x 1/4"	43 x 27.9 x .6 cm
40404	Large	22 x 20 x 1/2"	56 x 50 x 1.3 cm



Patient Axillary Positioner

AORN encourages the use of Axillary Rolls in place of sandbags or IV bags.

These extremely versatile positioners are good for many procedures.

40602	Small	10 x 4 1/4 x 2"	25.4 x 11.4 x 5 cm
40605	Large	14 x 4 x 2"	35.5 x 10.2 x 5 cm
40615	Extra Long	18 x 4 x 2"	46 x 10 x 5 cm



Heel Support

Place patient's heel in cup cut-out. Slit in back of cup secures the patient's Achilles tendon area.

40502	6 x 3 1/4" x 1 3/4"	15 x 8.3 x 4.5 cm
--------------	---------------------	-------------------



Positioning Strap

Welded seams and waterproof film combined with a soft patient-facing polymer surface make this a must have for the operating room. Radiolucent and MRI safe. Hook and loop attachments secure the strap to standard tables. The longer strap was specifically designed for larger patients. *This strap is not meant to restrain patient.*



Resurfacing Kits

Quick and easy resurfacing of minor cuts or abrasions on any Akton® polymer surface.

Watch our video on resurfacing your product.
<https://www.actionproducts.com/resources/downloads/videos.html>

KIT002	4 x 4"	10 x 10 cm
KIT003	36 x 4"	91.4 x 10 cm

40713	57 x 4 x 1/4"	144.8 x 10.2 x .6 cm	Adjustable range: 28 x 43"	71.12 x 109.22 cm
40713L	90 x 4 x 1/4"	228.6 x 10.2 x .6 cm	Adjustable range: 59 1/2 x 79 1/2"	151.13 x 201.93 cm



From Pediatrics to Bariatrics, we have your safe patient positioning needs covered.

Action Products also manufactures and sells a full line of Patient positioners made with our proprietary Akton® polymer.



Action Products, Inc.

954 Sweeney Drive
Hagerstown, MD 21740, USA

Customer Service:

Monday - Friday: 8:30am to 5:00pm EST

P 800.228.7763 or 301.797.1414

F 877.733.2073 or 301.733.2073

service@actionproducts.com

www.actionproducts.com



www.actionproducts.com/blogs



Product Warranty

Action Products warrants its standard line of Action® Operating Room Pads and Positioners for two years for defects in materials and workmanship at the time of purchase with the exception of products that are obsolete, misused, abused, altered, damaged, specialty or custom manufactured. Accessory and disposable items are warranted to be free of defective materials and workmanship at the time of initial purchase only.

General Terms

(Please note contracted terms may differ slightly) New accounts subject to credit approval. Standard credit terms are Net 30. We accept Visa®, MasterCard®, and American Express®. Standard domestic shipments are made by our preferred carrier, UPS, FOB Action Products, prepaid and added. Handling and drop-ship fees may apply. Specifications, prices, terms and availability subject to change without notice.

Product Returns

All product returns must be pre-authorized by a Customer Service Representative who will issue a unique RMA (Return Merchandise Authorization) number. Customer must supply order number and reason for return information. Returns are subject to a 20% restocking fee unless related to a manufacturing defect, product recall, product withdrawal, product correction or Customer Service error. Action Products will not accept order cancellations or returns of custom manufactured items.

The RMA number must accompany the returned item(s) and be clearly marked on the outside of the return shipping carton.

ACTION, AKTON and Making A Material Difference are registered trademarks of Action Products, Inc. USA

Position of model is intended to display pads and may not be surgically correct in every case. All product sizes are subject to change slightly.

© 2023 Action Products, Inc.