



Making a Material Difference®

Action Products, Inc.
www.actionproducts.com

INTERNATIONAL CATALOG

Patient Positioning from Pediatrics to Bariatrics

Clinically proven Action O.R. Support Surfaces and Positioners

Leg, Ankle &
Foot Positioning
p. 11




Overlays & Disposable
Overlay Covers
p. 3

Head & Neck
Positioners
p. 4

Supine Positioning
p. 18

Arm, Shoulder &
Elbow Positioners
p. 6

Chest, Torso & Hip Positioners
p. 8

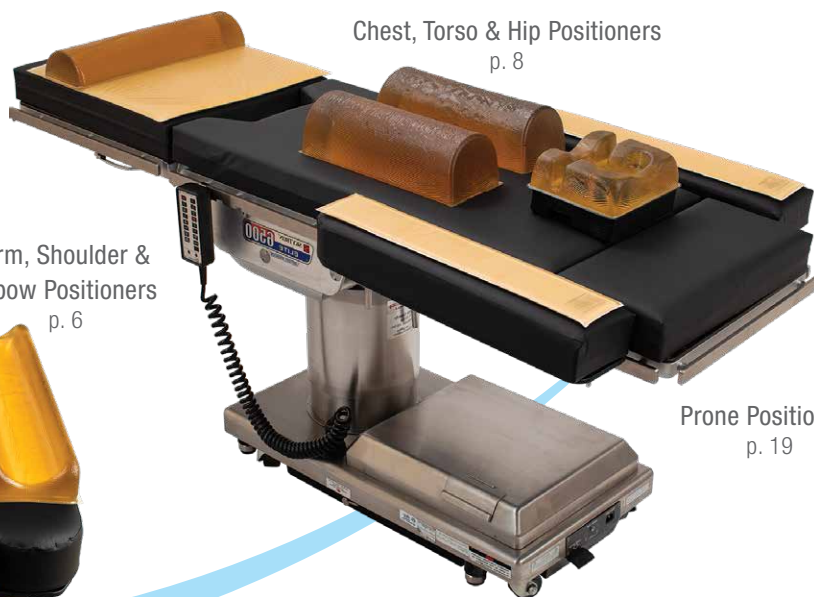

Prone Positioning
p. 19

Pediatric & Neo-Natal
Positioners
p. 16



Akton® Polymer Benefits p. 2



AMERICAN MADE
SINCE 1970



Making a Material Difference®

Advantages of AKTON viscoelastic polymer


- Excellent shear reducing qualities
- Soft, flexible film
- Enhances patient stability
- Evenly distributes heat and pressure
- Clinically proven to help prevent pressure injuries
- Reusable & Cost Effective
- Does not support bacterial growth
- Hypoallergenic
- Imaging compatible
- Biocompatible
- Hypo-allergenic
- Not made with natural rubber latex
- Free of phthalates, plasticizers & silicone
- Free of DEHP




Our reusable pads and positioners are impervious to fluids, will not support bacterial growth, do not leak or flow if punctured, and are easy to clean and maintain.



The shear reducing qualities of the products are enhanced by the combination of the soft, flexible outer film and the tissue equivalent polymer.

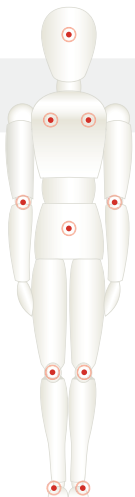

FACTORS CONTRIBUTING TO PRESSURE ULCERS:

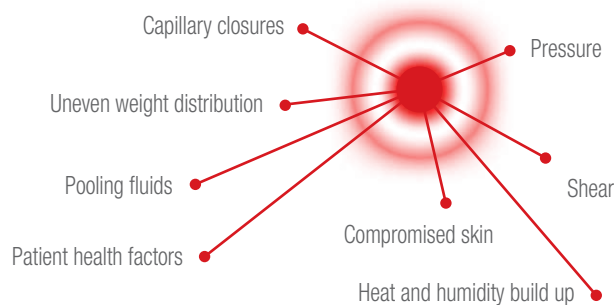

How our proprietary polymer reduces the risk...

The ultrasoft polymer combined with the soft, flexible surface of our pads allows the product to move with the patient, reducing negative effects of shear on fragile skin. In addition, the Akton® polymer distributes weight and dissipates heat; reducing the risk of high-pressure points and hot spot development that can lead to pressure injuries. These combined benefits enhance patient protection.

Look for authentic Akton® polymer

► If you have questions, email Customer Service at service@actionproducts.com




Overlays

The Action® O.R. Overlay is a versatile polymer pad used for a wide range of procedures. Its shear reducing qualities combined with the ability to distribute heat and pressure make it a valuable pressure management tool.

Used in longer procedures or for patients with fragile skin integrity, the overlay is an affordable upgrade for basic foam table pads that provides enhanced patient protection. The overlays are portable and can be exchanged between rooms as needed and are available in all standard OR table footprints in 1/2" thickness. Disposable covers are available for select overlay models.

40105	Small	20 L x 20 W x 1/2 H"	50.8 x 50.8 x 1.3 cm
40100	Medium	46 L x 20 W x 1/2 H"	117 x 50 x 1.3 cm
40109	Medium Large	63 L x 20 W x 1/2 H"	160 x 50 x 1.3 cm
40101	Full Length	72 L x 20 W x 1/2 H"	183 x 50 x 1.3 cm
40103	Perineal Cutout	46 L x 20 W x 1/2 H" w/ cutout	117 x 50.8 x 1.3 cm w/ cutout
40104	3-piece Set	See A, B and C	
40104A	Head	9 L x 20 W x 1/2 H"	22.9 x 50.8 x 1.3 cm
40104B	Torso w/ perineal cutout	39 1/2 L x 20 W x 1/2 H"	100.3 x 50.8 x 1.3 cm
40104C	Foot w/ perineal cutout	20 L x 20 W x 1/2 H"	50.8 x 50.8 x 1.3 cm

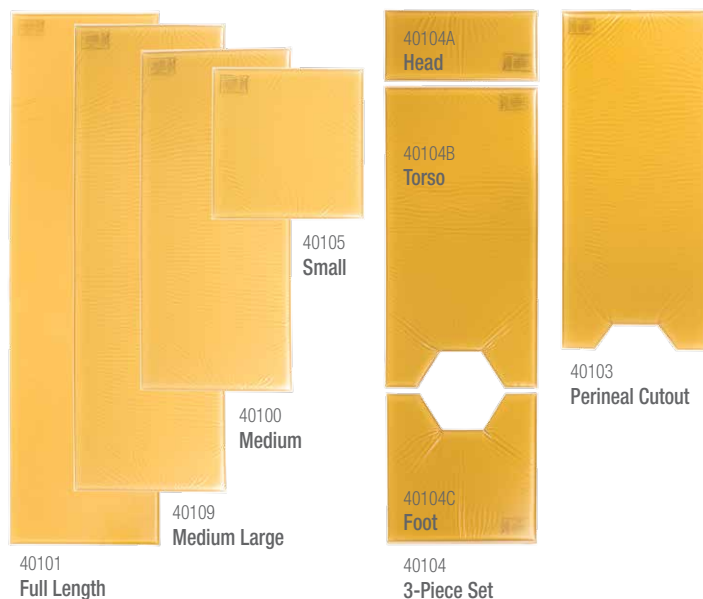

AORN 2025 Guidelines, Pressure Injury Prevention

Selecting and Using Support Surfaces, section 3.2

“Position patients on surfaces that reduce the potential for PI by redistributing pressure and reducing the risk of unintended patient movement.... Static surfaces such as viscoelastic polymer (i.e., gel) overlays are designed to reduce shearing and to support the patient's weight without becoming fully compressed under the weight of the patient's body.^{5,46,74} Gel overlays are radiolucent, latex-free, and reusable.”⁷⁴


Disposable Overlay Covers

Action's fluid-resistant, non-sterile Disposable Overlay Covers are made of a soft, flexible film and are available in two styles. They are designed to help reduce exposure to fluids and to help extend the life of the overlay.

The “pillowcase” style is closed at one end and fully encases the overlay like a sheath.

The “fitted sheet” style covers the top of the pad and is held in place by two pockets at each end of the cover. This cover has a 4" film drape that runs the length of the cover on both sides.

Fitted Sheet Style		
40800	Made for overlay model 40100- medium overlay	Case of 25
40801	Made for overlay model 40101- full length overlay	Case of 25
40803	Made for overlay model 40103- medium with cut out	Case of 25
40805	Made for overlay model 40105- small overlay	Case of 25
Pillowcase Style		
40102	Fits overlay models 40100 & 40101	Case of 25

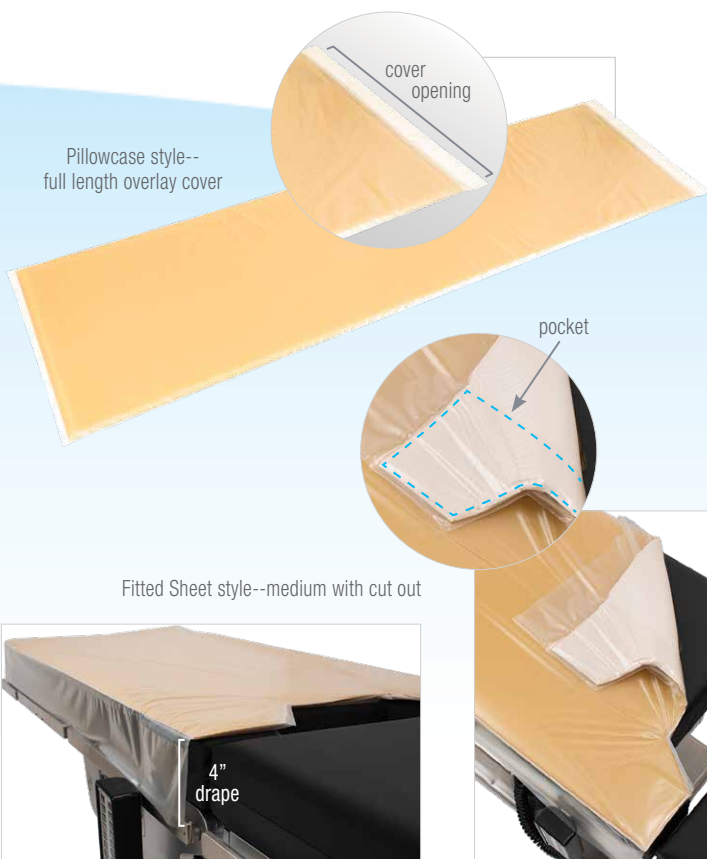

Donut Head Pads

Designed to protect and cradle the patient's head, these pads afford easy access for surgical area of the head and neck. It can be used for a variety of procedures including ENT, plastic, general and ophthalmic.

40201	Adult	8 OD x 3 ID x 1 ¾ H"	20.3 OD x 7.6 ID x 4.5 cm
40212	Adult, High	8 OD x 3 ID x 3 H"	20.3 OD x 7.6 ID x 7.6 cm
40218	Adult, Wide	8 OD x 4 ½ ID x 1 ¾ H"	20.3 OD x 11.4 ID x 4.5 cm
40202	Pediatric	5 ½ OD x 2 ¼ ID x 1 ¼ H"	14 OD x 5.7 ID x 3.2 cm
40219	Pediatric/ Neo-Natal	4 ½ OD x 2 ID x 1 H"	11.4 OD x 5.1 ID x 2.5 cm
40203	Neo-Natal	3 ¼ OD x 1 ½ ID x ¾ H"	8.3 OD x 3.8 ID x 2 cm



Donut Head Pads with Center Dish

40201D	Adult w/ Center Dish	8 OD x 3 ID x 1 ¾ H"	20.3 OD x 7.6 ID x 4.5 cm
40202D	Pediatric w/ Center Dish	5 ½ OD x 2 ¼ ID x 1 ¼ H"	14 OD x 5.7 ID x 3.2 cm
40203D	Neo-Natal w/ Center Dish	3 ¼ OD x 1 ½ ID x ¾ H"	8.3 OD x 3.8 ID x 2 cm



Supine Head Pad

The Supine Head Pad features a dual layer foam and polymer construction. The supportive foam base has a 6" circular cutout positioned beneath the soft, flexible polymer pad.

The foam cutout allows the 1" moderately dished polymer pad to immerse the occiput into the pad, creating a cradling effect that increases surface contact, lowers pressure, and gently positions the head. The entire head pad is encased in a fluid-resistant, flexible film cover to protect the porous foam base from contaminants and should not be removed.



40200	10 L x 9 W x 2 ¼ H	25.4 x 22.9 x 5.7 cm
--------------	--------------------	----------------------

Contoured Head Pad

The contoured design provides stability of the patient's head while in the supine position.



40213	7 ½ L x 7 ½ W x 3" H (1 ½" contoured)	19.4 x 19.4 x 7.6 cm
--------------	---------------------------------------	----------------------



Head Pad with Center Dish

This moderately dished head pad designed for lateral position reduces pressure on the side of the patient's face and ear.

40207	Adult	9Lx10Wx1"H	22.9x25.4x2.5 cm
40221	Pediatric	6Lx6Wx1"H	15.2x15.2x2.5 cm



Adult




Pediatric

Forehead Protector

The Forehead Protector is semi-circular with a flat bottom. It is used to protect the infant, child or adult patient's forehead from pressure exerted by oxygen tubes.




40208	2 ½ L x 1 ½ W x ¾ H	6.4 x 4.1 x 1.6 cm
--------------	---------------------	--------------------

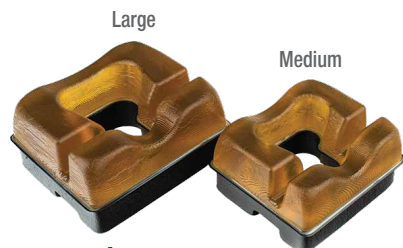

Prone Headrests


Our popular Prone Headrest is designed to support the patient's head in a straight, neutral position during prone procedures. Sculpted central cutout is designed to avoid pressure around eyes, nose, and face. Soft, Akton® polymer surface helps reduce shear at the forehead and chin support area. Channels present in top and base of pad provide tube access for the anesthesiologist.

40216M	Medium	8 ¾ L x 7 ¾ W x 4 ¾ H	22.2 x 19.7 x 12 cm
40216L	Large	10 ½ L x 9 W x 5 ¼ H	27 x 22.9 x 14 cm

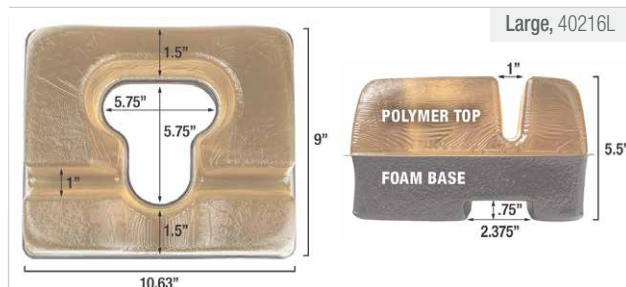

Ophthalmic Cradle Headrest

The unique contours of this supine head pad offer broad support for the neck of the patient while positioning the head in an oval center dish to reduce pressure and improve head stability during use. Suitable for most supine procedures, the head pad was originally designed for longer ophthalmic surgeries and features two raised areas on each side of the pad to offer soft support for a surgeon's hands. Well-suited option for patients requiring additional neck support.

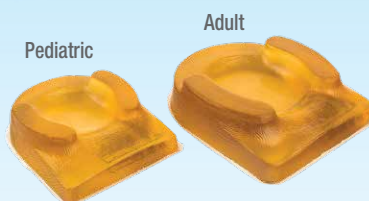

40217	Adult	10 ½ L x 9 W x 2 ⅝ H	26.7 x 22.9 x 6.7 cm
40217P	Pediatric	8 ½ L x 7 ½ W x 2 ¼ H	21.6 x 19 x 5.7 cm

Ophthalmic Headrests

This Ophthalmic Headrest features a deep contoured cutout that cradles the head and helps enhance stability. Its robust, soft sides are designed to lessen surgeon fatigue and offer wrist and hand support during ophthalmic, facial, or oral procedures.



40209	Large	10 L x 12 W x 3 ⅞ H	25.4 x 30.5 x 9.8 cm
40210	Medium	10 L x 12 W x 2 ⅞ H	25.4 x 30.5 x 7.3 cm

Horseshoe Head Pads

The horseshoe configuration provides a clear air passageway and aids the anesthesiologist while patients are in a lateral and prone position, head to side.

40204	Adult	8 OD x 3 ID x 1 ¾ H	20.3 OD x 7.6 ID x 4.5 cm
40211	Adult, High	8 OD x 3 ID x 3 H	20.3 OD x 7.6 ID x 7.6 cm
40205	Pediatric	5 ½ OD x 2 ¼ ID x 1 ¼ H	14 OD x 5.7 ID x 3.2 cm
40220	Pediatric/ Neo-Natal	4 ½ OD x 2 ID x 1 H	11.4 OD x 5.1 ID x 2.5 cm
40206	Neo-Natal	3 ¼ OD x 1 ½ ID x ¾ H	8.3 OD x 3.8 ID x 2 cm






Hand / Wrist Support Pad

Reusable, polymer pad offers stable and comfortable hand and wrist positioning and replaces rolled up towels. Perfect for Carpal Tunnel, Arterial Cannulation, Dialysis and IV or A-line insertion support.

40321 11 L x 3 ½ W x 2 ¼ H" 27.94 x 8.89 x 5.7 cm



Shown with Stabilizer Strap

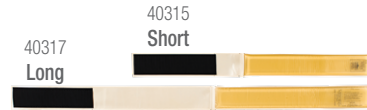

Arm / Hand Stabilizer Strap

This reusable strap features a soft polymer section to protect patient's skin while stabilizing arm/ hand. Strap is secured with attached hook and loop.

40315	Short	25 L x 2 W x ¼ H" 12 ½ L x 2 W x ¼ H" (polymer area only)	63.5 x 5.08 x .64 cm
40317	Long	36 L x 2 W x ¼ H" 15 ½ L x 2 W x ¼ H" (polymer area only)	91.44 x 5.08 x .64 cm

Product shown in use on Hand/Wrist Support, Armboard Pads, and Contoured Armboard Pads.



Shown with Stabilizer Strap



40307 Small 40310 Medium 40306 Large

Contoured Armboard Pads

The contoured armboard design cradles patient's arm while providing full-arm protection of the bicep, elbow, and forearm. Gently raised sides to help keep arm securely positioned.

40307	Small	16 L x 5 ¾ W x 2 ¼ H"	40.64 x 14.61 x 5.72 cm	3/4" Deep Middle Channel
40310	Medium	20 L x 5 ¾ W x 2 ¼ H"	50.8 x 14.6 x 5.72 cm	
40306	Large	24 L x 5 ¾ W x 2 ¼ H"	60.96 x 14.6 x 5.72 cm	



Shown with Stabilizer Strap



40301 Long 40308 Wide Long 40300 Short 40309 Wide Short

Armboard Pads

For use on standard armboard table extensions to protect the patient's bicep, elbow and forearm areas.

40300	Short	15 L x 4 W x ¼ H"	38.1 x 10.16 x .64 cm
40301	Long	24 L x 4 W x ¾ H"	60.96 x 10.16 x .95 cm
40309	Wide Short	15 L x 6 ½ W x ¼ H"	38.1 x 16.51 x .64 cm
40308	Wide Long	24 L x 6 ½ W x ¾ H"	60.96 x 16.51 x .95 cm

Ulnar / Forearm Protector

This pad helps to protect the patient's ulnar nerve and full forearm from damage due to pressure or shear.

Designed to be wrapped around the arm and secured with the attached hook & loop.

40303 18 L x 13 W x ¼ H" 45.72 x 33.02 x .64 cm



Hook and loop attachment

Toboggan Pad

This pad is for use with the arm/ leg guard shield (toboggan). Together they secure the arm or leg position on the table and protect the patient from being leaned on or pinched during the procedure.

The polymer pad cushions and helps protect the arm's ulnar and radial nerves.

The small film pocket slips on to the shield and pad is secured by hook & loop attachment (included).

The shield is only sold with the set.



Shown with Heel Support



Hook and loop attachment



40130	Toboggan Pad	16 ¼ L x 8 ¼ W x ¾ H"	41.27 x 20.96 x .95 cm
S40130*	Toboggan Pad Set	16 ¼ L x 8 ¼ W x ¾ H"	41.27 x 20.96 x .95 cm
	Shield size =	15 ½ L x 18 ¼ W"	39.37 x 46.4 cm


*Warranty covers pad only.

*Toboggan Armguard Shield is not available for sale or shipment outside of the United States and Canada.

Mayo Stand Cover

Use to prevent pressure and nerve damage when limb is positioned on the Mayo Stand during procedures.

5/8" polymer pad is sealed along three sides to a thick, durable film creating a pocket that slides easily on to the top of the Mayo stand.

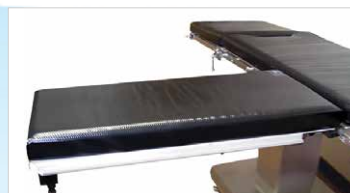

AKTON polymer with soft, flexible film surface



High-density Foam

40717 21 L x 18 ½ W x ¾ H" 53.3 x 47 x 1.6 cm

“The Core Curriculum for the RN First Assistant, 5th Edition states that, *during the Shoulder Replacement, the operative arm will be draped free for intraoperative manipulation; mechanical devices may be used to assist with manipulation or a padded Mayo stand can be used to support the arm.*

“Additionally, AORN Guideline 3.15 states: *Prevent the patient's body from contacting metal portions of the OR bed and other hard surfaces (eg, the Mayo Stand).* [Recommendation]

Hourglass Hand Table Pad

Reusable 1/2" polymer overlay designed for surgical procedures involving the forearm, elbow, or hand utilizing the hourglass hand table. Good for lengthy hand microsurgeries requiring additional pressure and shear relieving support. Fits most standard hourglass hand tables.



40330 32 L x 24 W x ½ H" 81.28 x 60.96 x 1.27 cm
16"W middle section 40.6 cm middle section

Hand Table Composite Pad

Designed for surgical procedures involving the forearm, elbow, or hand. Good for lengthy hand microsurgeries requiring additional pressure and shear relieving support. Pad combines an AKTON polymer patient surface atop a high-density foam base featuring fully-welded, non-stitch construction. This pad is designed/constructed with the same pressure management benefits as the Action Products' "Original" O.R. table pads.

40320 34 L x 16 W x 2 ¾ H" 86.36 x 40.64 x 6.03 cm



Ulnar Nerve Protector

This pad helps to protect the patient's ulnar nerve while providing greater visibility and access for IV or A-line for anesthesiologist. Pad is placed around the elbow and secured with attached hook & loop. Maximum pad circumference is approximately 15.75" on the small and 19" on the large.

40304	Small	15 ¾ L x 6 W x ½ H"	40 x 15.24 x 1.27 cm
40304L	Large	19 L x 8 ½ W x ½ H"	48.26 x 21.6 x 1.27 cm
40304P	Pediatric	10 L x 3 W x ¼ H"	25.4 x 7.6 x .6 cm



Pediatric



Chest Rolls

This line of flat bottom chest rolls is designed to provide stable support in the prone position. They are available in a wide range of heights and lengths to accommodate virtually any patient, adult to neonate.

Adult		
40616	12 ½ L x 4 ¾ W x 3 ¾" H	31.75 x 11.11 x 8.57 cm
40617	16 ½ L x 4 ¾ W x 3 ¾" H	41.91 x 11.11 x 8.57 cm
40618	20 ½ L x 4 ¾ W x 3 ¾" H	52.07 x 11.11 x 8.57 cm
40619	12 L x 6 W x 5" H	30.48 x 15.24 x 12.7 cm
40620	16 L x 6 W x 5" H	40.64 x 15.24 x 12.7 cm
40621	20 L x 6 W x 5" H	50.8 x 15.24 x 12.7 cm

Pediatric		
40639	6 L x 2 W x 2" H	15.24 x 5.08 x 5.08 cm
40640	8 L x 2 W x 2" H	20.32 x 5.08 x 5.08 cm
40641	10 L x 2 W x 2" H	25.4 x 5.08 x 5.08 cm
40642	12 L x 2 W x 2" H	30.48 x 5.08 x 5.08 cm
40643	6 L x 3 W x 3" H	15.24 x 7.62 x 7.62cm
40644	8 L x 3 W x 3" H	20.32 x 7.62 x 7.62cm
40645	10 L x 3 W x 3" H	25.4 x 7.62 x 7.62cm
40646	12 L x 3 W x 3" H	30.48 x 7.62 x 7.62cm



Disposable Chest Roll Covers

Action's fluid-resistant, non-sterile Disposable Chest Roll Covers are made of a soft, flexible film designed to extend the life of the chest roll and reduce exposure to fluids during use.

The slip-on covers mimic size and shape of the chest roll and have an overlapping film flap along the flat bottom.

Purchase in cases of 25.

Fits Action® Chest Roll		
40616DC	40616	Case of 25
40617DC	40617	Case of 25
40618DC	40618	Case of 25
40619DC	40619	Case of 25
40620DC	40620	Case of 25
40621DC	40621	Case of 25



Axillary Roll Cover

Designed to upgrade basic foam positioning rolls with reusable polymer wrap-around covers. Covers are adjustable and are secured by hook and loop attachment.

40601 17 L x 12 W x ½" H 43.18 x 30.48 x 1.27 cm



Chest Roll Cover

Designed to upgrade basic foam chest rolls with reusable polymer wrap-around covers. Covers are adjustable and are secured by hook and loop attachment.


40600 18 L x 18 W x ½" H 45.72 x 45.72 x 1.27 cm



Lateral Positioner

Designed to protect the patient while in the lateral position. Accommodates both shoulder and upper arm. Provides pressure relief and nerve protection. Includes a lightweight foam base with AKTON polymer pad overlay.

40610 28 L x 19 1/2 W x 6 - 1" H 71.12 x 49.53 x 15.24 - 2.54 cm

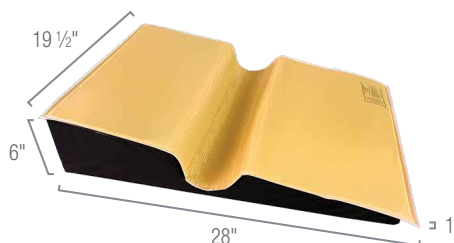


Patient Positioner

AORN encourages the use of Axillary Rolls in place of sandbags or IV bags. These extremely versatile positioners are good for many procedures.

40602 Small 10 L x 4 1/4 W x 1 1/8" H 25.4 x 10.8 x 4.76 cm

40605 Long 14 L x 4 W x 2" H 35.56 x 10.16 x 5.08 cm

40615 Extra Long 18 L x 4 W x 2" H 45.72 x 10.16 x 5.08 cm

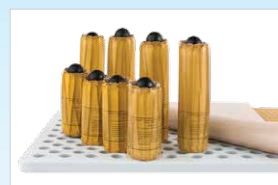

Peg Board Pad Set

The Action® Peg Board Pad set complements most standard peg board positioning systems. This multi-piece set is designed to protect patient skin and soft tissue from damaging pressure and shear forces while in the high-pressure lateral position.

The medium-length overlay (Item #40100) lies atop the peg board (beneath the patient) while the adjustable peg post pads of varying lengths easily wrap around the pegs and are secured by hook and loop attachment. Overlay and peg pads are available as a set, or as individual units.

Set (#40719) includes:


- 1 Medium Overlay Pad, #40100
- 4 Short Peg Post Pads, #40719A
- 4 Tall Peg Post Pads, #40719B

Peg board and peg posts not included

Graduated Thoracic Positioner

The graduated, or sloping, configuration of this pad encourages full chest extension. It is designed for cardiovascular or thoracic procedures, which mandate full chest hyperextension.

40604 14 L x 9 W x 2 - 1 1/2" H 35.56 x 22.86 x 5.08 - 3.81 cm

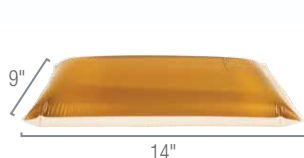


Chest Positioning Set

These positioners have a contoured design patterned to uniformly conform to a patient's shoulder to chest area, with shaped edges for the axilla. They provide upper body positioning while alleviating pressure on the thorax area.

40625 Left 28 - 23 3/8 L x 8 x 3 x 2 1/2" H 71.12 x 58.42 x 20.32 x 7.62 x 6.35 cm

40626 Right 28 - 23 3/8 L x 8 x 3 x 2 1/2" H 71.12 x 58.42 x 20.32 x 7.62 x 6.35 cm

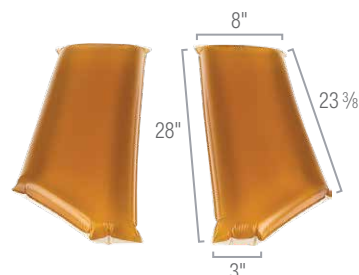

Positioning Strap

Welded seams and fluid-resistant film combined with a soft patient-facing polymer surface make this a must have for the operating room.

Radiolucent and MRI safe.

Hook and loop attachments secure the strap to standard tables. The longer strap was specifically designed for larger patients.

This strap is not meant to restrain patient.



40713	57 L x 4 W x ¼" H	144.78 x 10.16 x .64 cm
	AKTON polymer surface measures 98 x 8.2 x 0.6cm	
	Adjustable range: 28 x 43" 71.12 x 109.22 cm	
40713L	90 L x 4 W x ¼" H	228.6 x 10.16 x .64 cm
	AKTON polymer surface measures 190.5 x 8.2 x 0.6cm	
	Adjustable range: 59½ x 79½" 151.13 x 201.93 cm	



Orthopedic Hip Positioner

Designed to complement the Anterior Supine Approach to THR, this unique hip positioner was engineered to provide consistent, optimal site access and patient stability while protecting the patients skin integrity.

41740	8 ½ L x 16 W x 2.63" H	21.59 x 40.64 x 6.68 cm
--------------	------------------------	-------------------------



Shoulder Brace, 2-Piece Set

Reusable 1" thick polymer pads designed for basic shoulder braces. Adjustable pads wrap around braces and are secured by hook and loop attachments. Sold as a pair.

40302	5 L x 6 W x 1" H	12.7 x 15.24 x 2.54 cm
--------------	------------------	------------------------

Vac-Pac® Covers




This pad is designed to line the interior of the Olympic Vac-Pac® (bean bag) positioner to help mitigate patient pressure and shear while in the lateral position. It is applicable for orthopedic, thoracic or urology procedures that dictate the patient be stable in the lateral position.

Vac-Pac® is a registered trademark of Natus Medical Incorporated

40612	Vac-Pac® Cover Medium Size 30	37 L x 28 W x ½" H	93.98 x 71.12 x 1.27 cm
40612-31	Vac-Pac® Cover Large Size 31	35 ½ L x 38 ¾ W x ½" H	90.17 x 98.43 x 1.27 cm
40613	Vac-Pac® Cover T-Shape Size 35	15 L x 34 W x ½" H (bottom) 9 L x 23 W x ½" H (top)	38.1 x 86.36 x 1.27 cm (bottom) 22.86 x 58.42 x 1.27 cm (top)

Breast Biopsy Table Pad

Improve patient comfort with this soft, polymer pad that helps protect skin and bony prominences while in the prone position during breast biopsy procedures.

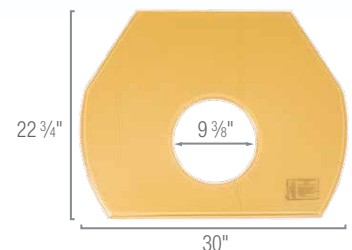

40715	30 L x 22 ¾ W x ¾" H with 9 ¾" D cutout	76.2 x 57.79 x .95 cm with 23.81 cm D cutout
--------------	---	--

Resurfacing Kits

Quick and easy resurfacing of minor cuts or abrasions on any Akton® polymer surface.

KIT002	4 x 4"	10.16 x 10.16 cm
KIT003	36 x 4"	91.44 x 10.16 cm





Heel Support

Reusable polymer heel support designed to float heels above table pad surface while in supine position.

Soft heel support features a v-shaped notch to position and accommodate the Achilles tendon and sculpted heel cup to cradle and securely position the foot in a neutral position.

40502 6 L x 3 1/4 W x 1 3/4" 15 x 8.3 x 4.5 cm



Heel Protector Pad

This "clamshell" style pad provides full coverage ankle and heel protection during lengthy surgeries. Interior surface of the pad is permanently labeled to ensure proper foot placement. Easy to use and secure with adjustable hook and loop attachment.

Accommodates approximately, a 15" ankle circumference. Use alone in supine and lateral position, or combine with Ankle Strap Stirrup and Candy Cane pads for lithotomy position.

40401 12 1/2 L x 7 3/4 W x 1/2 H" 32 x 20 x 1.2 cm



Knee Crutch Pads

Three sizes and thicknesses available to line knee crutches and help reduce pressure and shear in the popliteal space.

SMALL: Thinnest pad, best for lightweight support for urology, OB/ GYN, and pediatric procedures.

MEDIUM: 85% larger than the Small, this 1/4" pad offers midweight support for urology and OB/GYN procedures.

LARGE: Heavyweight 1/2" thick, 22 x 20" pad best suited for large knee crutches.



40402	Small	11 L x 9 W x 3/8 H"	27.9 x 22.9 x .95 cm
40403	Medium	17 L x 11 W x 1/4 H"	43 x 27.9 x .6 cm
40404	Large	22 L x 20 W x 1/2 H"	56 x 50 x 1.2 cm



Dome-Shaped Positioners

Shown with Heel Support and Overlay

Semi-circular design with flat bottom configuration provides pressure protection under the knees while the patient is in a supine position. Helps to protect patient against foot drop for patients in prone position.

Can be used under the metatarsal or under the foot, against the board, or in the semi-Fowler position.

40603	14 L x 6 W x 3 H"	36 x 15 x 7.6 cm
40603L	20 L x 6 W x 3 H"	50.8 x 15 x 7.6 cm

Temporarily unavailable

Positioning Strap

Welded seams and fluid-resistant film combined with a soft patient-facing polymer surface make this a must have for the operating room.


Radiolucent and MRI safe.

Hook and loop attachments secure the strap to standard tables. The longer strap was specifically designed for larger patients.

This strap is not meant to restrain patient.



40713	57 L x 4 W x ¼" H	144.78 x 10.16 x .64 cm
	AKTON polymer surface measures 98 x 8.2 x 0.6cm	
	Adjustable range: 28 x 43" 71.12 x 109.22 cm	
40713L	90 L x 4 W x ¼" H	228.6 x 10.16 x .64 cm
	AKTON polymer surface measures 190.5 x 8.2 x 0.6cm	
	Adjustable range: 59½ x 79½" 151.13 x 201.93 cm	

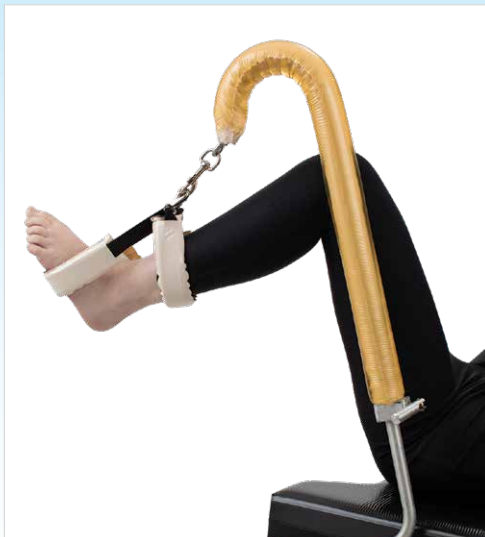

Candy Cane Stirrup Pad Sets

Use in lithotomy position during gynecological and urological procedures; protects peroneal nerve from pressure and shear injuries that can result from leg contacting metal stirrups. Soft, durable and reusable product eliminates need for tedious wrapping and rewrapping of the leg holder posts.

Candy Cane Pad Set



Recommend using Action® Ankle Stirrup Pad Set and full coverage Heel/ Ankle Protectors for enhanced patient protection.

**Lithotomy holders not included*

40501	Candy Cane (Set of 2)	28 L x 2 OD x 1 ID" 71 x 5 x 2.5 cm
--------------	-----------------------	-------------------------------------

Ankle Strap Stirrup Pad Sets

Two adjustable polymer straps are combined to support and protect the patient's leg and foot while leg is suspended in lithotomy position.

Ankle Strap Stirrup Pads

Large



Small



40510	Ankle Strap Stirrup Pads, Small (Set of 2)	12 ½ L x 1 ¾ W x ¾ H" 31.7 x 4.7 x .95 cm
40511	Ankle Strap Stirrup Pads, Large (Set of 2)	17 ¾ L x 1 ¾ W x ¾ H" 44.1 x 4.7 x .95 cm

Heel Comfort Cup Set

Lightweight polymer heel cups for comfort and support protects heels from pressure and shear injury. Can be used in conjunction with anti-embolism stockings. Appropriate for recovery and radiology. Comfortable all-polymer fastening strap attaches with simple hook and loop. One size fits all.

40520	(Set of 2)	4 ¼ L x 3 ¾ W x ½ H" 10.8 x 9.2 x 1.2 cm
--------------	------------	--



Toboggan Pad

This pad is for use with the arm/ leg guard shield (toboggan). Together they secure the arm or leg position on the table and protect the patient from being leaned on or pinched during the procedure.

The polymer pad cushions and helps protect the arm's ulnar and radial nerves.

The small film pocket slips on to the shield and pad is secured by hook & loop attachment (included).

The shield is only sold with the set.



Shown with Heel Support



Hook and loop attachment



40130	Toboggan Pad	16 1/4 L x 8 1/4 W x 3/8 H"	41.27 x 20.96 x .95 cm
S40130*	Toboggan Pad Set	16 1/4 L x 8 1/4 W x 3/8 H"	41.27 x 20.96 x .95 cm
	Shield size =	15 1/2 L x 18 1/4 W"	39.37 x 46.4 cm

*Warranty covers pad only.

*Toboggan Armguard Shield is not available for sale or shipment outside of the United States and Canada.



One Piece



Two Piece Set



Frog Leg Positioner

Designed for saphenous vein harvesting; the contouring and construction of the pad enhance site access while providing stable support.


Akton® polymer overlay protects from pressure and shear injury and stable foam base (visco foam atop high-density foam) is encased in 4-way stretch water resistant material that is fully welded to prevent fluid ingress.

Reusable product eliminates need for wasteful use of stacked blankets or towels and sandbags.


Replacement polymer overlay is available, Item 40718LPA. Overlay features hook and loop to attach it to the fabric-covered base.

Two-piece design allows for independent leg positioning.

40718LP	Frog Leg Positioner	32 L x 21 W x 5 H" includes 1/4" AKTON polymer	81.3 x 53.3 x 12.7 cm includes 0.6 cm AKTON polymer
40718LPA	Frog Leg Polymer Replacement	35 L x 21 W x 1/4 H"	88.9 x 53.3 x .6 cm
40718LP-2	Frog Leg Positioner (2-Piece Positioner)	32 L x 10 1/2 W x 5 H" includes 1/4" AKTON polymer	81.3 x 26.6 x 12.7 cm includes 0.6 cm AKTON polymer

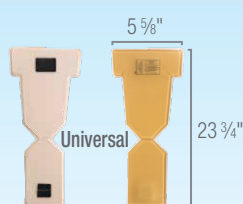

Universal

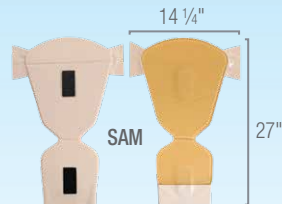

SAM

Boot Stirrup Pad Set

Soft, reusable 1/4" thick polymer lithotomy boot liner designed to protect patient's lower leg and heel area by reducing pressure and shear.



Shown with SAM Boot

For use with universal stirrups and includes hook and loop attachments to secure to hardware. Used in lithotomy position for gynecological and urological cases.

40500	Universal (Set of 2)	23 3/4 L x 5 5/8 W x 1/4 H"	60.3 x 13.7 x .6 cm
40503	SAM (Set of 2)	27 L x 14 1/4 W x 1/4 H"	68.5 x 36 x .6 cm

Fracture Table Foot Pad

For use on leg-prep stand and fracture table foot holders, for protection from pressure and shear injury.



40400	12 L x 4 W x 1/4 H"	30.5 x 10 x .6 cm
--------------	---------------------	-------------------



Our reusable products are engineered to help preserve skin integrity and provide comfort and stability for highly specialized surgeries.

Shown left: Wilson Frame Pad Set with Action® Prone Headrest, Composite Contoured Armboard, Small Overlay, and Chest Roll



Wilson Frame Pad Set

Designed to cover the Wilson Frame and protect patients from shear and pressure injury. These 1/2" thick pads, and robust 1" thick pads easily slip over frame and should be secured with Hook and Loop attachment.

For orthopedic and neuro procedures when Wilson Frame is used. Two-piece set.

40701	1/2" Thick	29 L x 10 W x 1/2 H"	74 x 25.4 x 1.3 cm
40702	1" Thick	29 L x 10 W x 1 H"	74 x 25.4 x 2.5 cm



Hook and loop attachment




Kambin Spinal Frame Pad Set

Designed to fit the original Kambin Spinal Frame.

This two-piece pad set fits over frame's patient supports. The 1/2" thick pads are secured to the frame by included hook & loop attachments positioned along the bottom of each pad as well as on tabs (2 each per pad) that fold over and attach to the outer edges of the frame.

40712	33 L x 13 3/4 W x 1/2 H"	13 L x 5.41 x 1.3 cm
--------------	--------------------------	----------------------

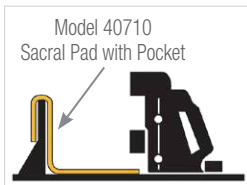

McGuire Hip/Sacral Pad

This 1/2" polymer pad protects the ilium, sacrum and area over the greater trochanter from shear and pressure damage while patient is in lateral position on the frame.

A pocket is created by folding the curved end of the pad over itself and securing attached hook and loop. To place on frame, slide pocket over the back plate, allowing rest of pad to cover the interior of the back plate and base of the frame.

40710	26 L x 10 W x 1/2 H"	66 x 22.5 x 1.3cm
--------------	----------------------	-------------------



Montreal Frame Pad

A positioning pad for the main platform of the Montreal Lateral Positioning Device. This pad is used in a lateral position to mitigate pressure and shear related to total hip procedures.

40110	18 1/4 L x 12 W x 1/2 H"	46.3 x 30.5 x 1.3 cm
--------------	--------------------------	----------------------

Fracture Table Hip Pad

Designed for use on the fracture table plate, the primary function of this pad is to reduce pressure and shear at the sacral area.



Fracture Table Hip Pad with Post Hole

This oval shaped hip pad features an internal 1.75" hole to accommodate the perineal post.

Pad is designed to protect the patient's upper leg and scrotum areas from skin injury due to pressure or shear forces during use.



Fracture Table Post Pad

This 1/2" thick wrap-around pad reduces pressure and shear on the patient's scrotum and upper leg areas.

It also eliminates the tedious chore of wrapping and rewrapping the fracture table post with various makeshift items such as gauze and towels.



40704
Post Pad

40703	Hip Pad	8 L x 5 1/2 W x 1/2 H"	20 x 14 x 1.3 cm
40704	Post Pad	7 L x 5 W x 1/2 H"	18 x 12.7 x 1.3 cm
40707	Hip Pad w/ Post Hole	3 1/4 x 3 3/8 x 1 1/4 OD"; 2 1/8 x 2 3/8 x 1 1/8 ID" 1 3/4" diameter hole to accommodate perineal posts	4.5 cm diameter hole to accommodate perineal posts



Andrews Frame™ Iliac Strap

Two-piece adjustable strap system featuring durable nylon fabric belt with removable polymer pad insert. Removable polymer pad reduces pressure and shear forces on iliac crest and upper thigh areas.

Polymer pad covers up to 24" across the hips. Strap is secured with adjustable hook and loop closure.

40706 Iliac Strap 24 L x 6 W x 3/8 H" 61 x 15 x .95 cm

Andrews Frame™ Knee Pad

Knee Pad designed for use with the Andrews Frame® on the knee platform to alleviate pressure and shear during use.

40705 Knee Pad 23 L x 16 1/4 W x 1 H" 56 x 40.6 x 2.5 cm

3/8" AKTON polymer pad



Action® Comfort Stool Cover

Action® Comfort Stool Cover - 3/8" thick Akton® polymer pad is bonded to a 4" nylon oxford band with drawstring closure.

40900 13 D x 3/8" H 33.02 x .95 cm

40901 14 D x 3/8" H 35.56 x .95 cm

40902 15 D x 3/8" H 38.1 x .95 cm

► Over 80% of the top U.S. pediatric hospitals use Action® products to protect their patients.

Donut Head Pads

Designed to protect and cradle the patient's head, these pads afford easy access for surgical area of the head and neck. It can be used for a variety of procedures including ENT, plastic, general and ophthalmic.

Available in traditional or dished style.

40202	Pediatric	5 ½ OD x 2 ¼ ID x 1 ¼ H"	14 OD x 5.7 ID x 3.2 cm
40219	Pediatric/ Neo-Natal	4 ½ OD x 2 ID x 1 H"	11.4 OD x 5.1 ID x 2.5 cm
40203	Neo-Natal	3 ¼ OD x 1 ½ ID x ¾ H"	8.3 OD x 3.8 ID x 2 cm

Teen / Adult sizes available



Pediatric Donut Head Pad



Donut Head Pads with Center Dish

40202D	Pediatric w/ Center Dish	5 ½ OD x 2 ½ ID x 1 ¼ H"	14 OD x 6.4 ID x 3.2 cm
40203D	Neo-Natal w/ Center Dish	3 ¼ OD x 1 ½ ID x ¾ H"	8.3 OD x 3.8 ID x 2 cm

Teen / Adult sizes available




Neo-Natal w/ Center Dish enlarged to show dished style.



Horseshoe Head Pads

The horseshoe configuration provides a clear air passageway and aids the anesthesiologist while patients are in a lateral and prone position, head to side.

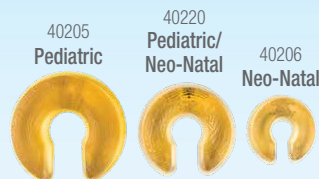

40205	Pediatric	5 ½ OD x 2 ¼ ID x 1 ¼ H"	14 OD x 5.7 ID x 3.2 cm
40220	Pediatric/ Neo-Natal	4 ½ OD x 2 ID x 1 H"	11.4 OD x 5.1 ID x 2.5 cm
40206	Neo-Natal	3 ¼ OD x 1 ½ ID x ¾ H"	8.3 OD x 3.8 ID x 2 cm

Teen / Adult sizes available

Head Pad with Center Dish

Designed for smaller patients, this moderately dished, square head pad may be used to reduce pressure on the occiput, neck, head or ears of the patient as needed.

40221	Pediatric	6 L x 6 W x 1 H"	15.2 x 15.2 x 2.5 cm
--------------	-----------	------------------	----------------------

Teen / Adult sizes available



Ophthalmic Cradle Headrest

The unique contours of this supine head pad offer broad support for the neck of the patient while positioning the head in an oval center dish to reduce pressure and improve head stability during use. Suitable for most supine procedures, the head pad was originally designed for longer ophthalmic surgeries and features two raised areas on each side of the pad to offer soft support for a surgeon's hands. Well-suited option for patients requiring additional neck support.

40217-P	Pediatric	8 ½ L x 7 ½ W x 2 ¼ H"	21.6 x 19 x 5.7 cm
----------------	-----------	------------------------	--------------------

Teen / Adult sizes available



Forehead Protector

The Forehead Protector is semi-circular with a flat bottom. It is used to protect the infant or child patient's forehead from pressure exerted by oxygen tubes.

40208	2 ½ L x 1 ¾ W x ¾ H"	6.4 x 4.1 x 1.6 cm
--------------	----------------------	--------------------





Ulnar Nerve Protector

This pad helps to protect the patient's ulnar nerve while providing greater visibility and access for IV or A-line for anesthesiologist. Pad is placed around the elbow and secured with attached hook & loop.

40304P Pediatric 10 L x 3 W x ¼ H" 25.4 x 7.6 x .6 cm

Teen / Adult sizes available



Universal Spot Protectors

This series of thin, round pads is designed to provide adaptive use padding options for those hard-to-protect areas unique to each patient. The shape is designed to enhance pressure redistribution and improve shear performance.

Thin profile, 1/4" thick polymer



40405	Mini	1 5/8" D x 1/4" H"	4.1 D x .6 cm	Pkg of 12
40406	Small	3 1/4" D x 1/4" H"	8.3 D x .6 cm	Pkg of 12
40407	Large	5 1/2" D x 1/4" H"	14 D x .6 cm	Pkg of 6

Bassinet Pad



NICU Pads

Provides pressure and shear management surface for pediatric and neo-natal patients.

40120L	Isolette Pad	26 L x 16 W x 1/2 H"	66 x 40.6 x 1.3 cm
40121	Isolette Pad	29 L x 19 W x 1/2 H"	74 x 48.2 x 1.3 cm
40120	Bassinet Pad	26 L x 12 W x 1/2 H"	66 x 30.5 x 1.3 cm



Pediatric Chest Rolls

Designed for smaller patients like neonates and toddlers, this line of flat bottom chest rolls provides stable support in the prone position. They are available in a wide range of heights and lengths to accommodate virtually any pediatric patient.

Teen / Adult sizes available

Pediatric/Neo-Natal Positioners

These positioners are used to position infants without rolling wash-cloths and other items.

40648	6 L x 1 1/2 W x 3/4 H"	15 x 3.8 x 2 cm
40649	6 L x 1 W x 1 H"	15 x 2.5 x 2.5 cm
40650	6 L x 1 1/2 W x 1 1/2 H"	15 x 3.8 x 3.8 cm

Teen / Adult sizes available

Patient Axillary Positioner

AORN encourages the use of Axillary Rolls in place of sandbags or IV bags. These extremely versatile positioners are good for many procedures.

40602	10 L x 4 1/4 W x 1 1/2 H"	25.4 x 10.8 x 4.76 cm
40605	14 L x 4 W x 2 H"	35.56 x 10.16 x 5.08 cm
40615	18 L x 4 W x 2 H"	45.72 x 10.16 x 5.08 cm



Pediatric Positioner Set


This versatile 4-piece set is designed to position and protect newborns, infants and toddlers during a variety of procedures. Placed flat and used alone, the channelled positioner helps stabilize patients while in supine. Used in combination with a chest roll under the pad at the thoracic area, the chest is extended for anterior chest procedures. Additionally, when used with the small chest roll under the channelled positioner at the thorax and large chest roll at the sacrum, access for posterior supine cases such as spinal tumor resection is enhanced.

Visit our website for more information on this set.



40614	Four piece pediatric positioning set Includes: 40205, 40642, 40646, 40647
40647	Rectangular Channel Pad 20 L x 9 W x 2 H" 50.8 x 22.9 x 5 cm

SUPINE Body Smart[®] Positioning Products

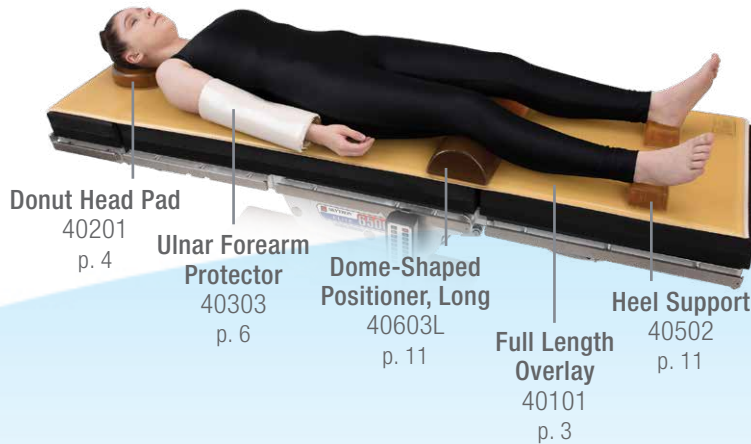

Head Positioners

- Donut Head Pad | p. 4
- Supine Head Pad | p. 4
- Head Pad with Center Dish | p. 4
- Contoured Head Pad | p. 4
- Ophthalmic Headrests | p. 5
- Ophthalmic Cradle Headrest | p. 5
- Forehead Protector | p. 4

Arm / Elbow / Hand Positioners

- Ulnar Forearm Protector | p. 6
- Ulnar Nerve Protector | p. 7
- Hand/Wrist Support | p. 6
- Contoured Armboards | p. 6
- Armboard Pads | p. 6
- Arm/Hand Stabilizer Straps | p. 6
- Toboggan Pad | p. 7
- Mayo Stand Cover | p. 7

Leg / Ankle / Heel Positioners

- Dome-Shaped Positioner | p. 11
- Positioning Strap | p. 12
- Heel Support | p. 11
- Heel Protector Pad | p. 11
- Heel Comfort Cup Set | p. 12
- Frog Leg Positioner | p. 13
- Axillary Roll Cover | p. 8

Support Surfaces

- Overlays | p. 3

LATERAL Body Smart[®] Positioning Products



Head Positioners

- Head Pad with Center Dish | p. 4

Arm / Elbow / Hand Positioners

- Contoured Armboards | p. 6
- Armboard Pads | p. 6
- Arm/Hand Stabilizer Straps | p. 6

Chest / Torso / Hip Positioners

- Patient Axillary Positioner | p. 9
- Peg Board Pad Set | p. 9
- Lateral Positioner | p. 9
- Positioning Strap | p. 12

Leg / Ankle / Heel Positioners


- Small Overlay | p. 3
- Heel Protector Pad | p. 11

Support Surfaces

- Overlays | p. 3

Positioning equipment, devices and support surfaces should be used in accordance with manufacturer's instructions for use (IFU). Hospital protocol will supersede above suggestions.

PRONE Body Smart[®] Positioning Products

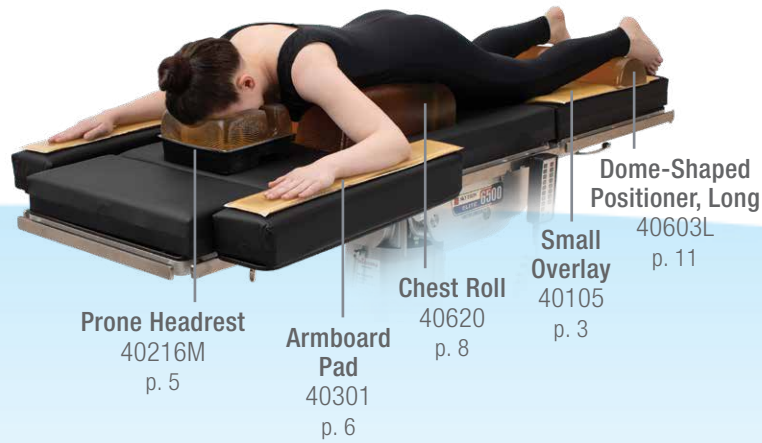

Head Positioners

- Prone Headrests | p. 5
- Horseshoe Head Pad | p. 5

Arm Positioners

- Armboard Pads | p. 6
- Arm/Hand Stabilizer Straps | p. 6
- Contoured Armboards | p. 6

Chest / Torso / Hip Positioners

- Chest Rolls | p. 8
- Chest Positioners | p. 9

Leg / Ankle / Heel Positioners

- Dome-Shaped Positioners | p. 11
- Axillary Roll Cover | p. 8


Specialty Frame Positioning

- Wilson Frame Pad Set | p. 14
- Kambin Spinal Frame Pad Set | p. 14

Support Surfaces

- Overlays | p. 3

LITHOTOMY Body Smart[®] Positioning Products

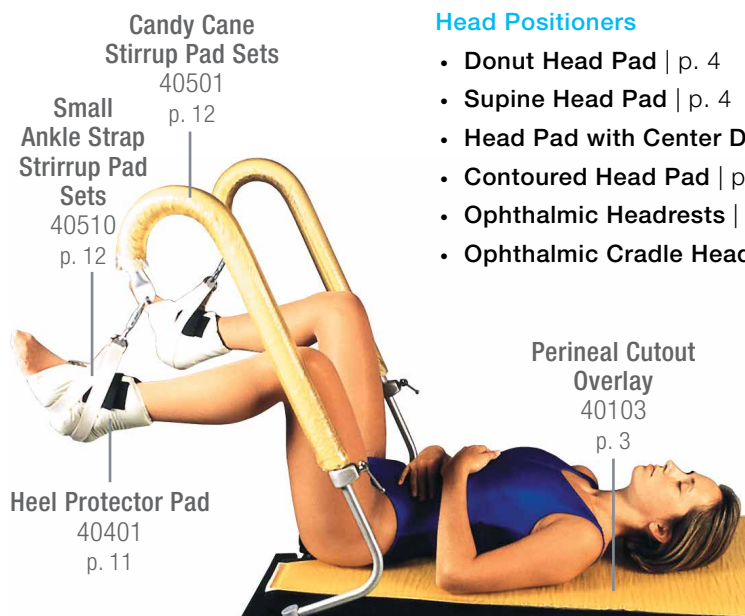

Head Positioners

- Donut Head Pad | p. 4
- Supine Head Pad | p. 4
- Head Pad with Center Dish | p. 4
- Contoured Head Pad | p. 4
- Ophthalmic Headrests | p. 5
- Ophthalmic Cradle Headrest | p. 5

Arm / Elbow / Hand Positioners

- Ulnar Forearm Protector | p. 6
- Ulnar Nerve Protector | p. 7
- Contoured Armboards | p. 6
- Armboard Pads | p. 6
- Arm/Hand Stabilizer Straps | p. 6

Leg / Ankle / Heel Positioners

- Candy Cane Stirrup Pad Sets | p. 12
- Ankle Strap Stirrup Pad Sets | p. 12
- Heel Protector Pad | p. 11
- Boot Stirrup Pad Set | p. 13
- SAM Boot Stirrup Pad Set | p. 13

Support Surfaces

- Overlays | p. 3



From Pediatrics to Bariatrics, we have your safe patient positioning needs covered.

Action Products manufactures and sells a full line of patient positioners
made with our proprietary Akton® polymer.

Action Products, Inc.

954 Sweeney Drive
Hagerstown, MD 21740, USA

Customer Service:

Monday - Friday: 8:30am to 5:00pm EST

P 800.228.7763 or 301.797.1414

F 877.733.2073 or 301.733.2073

service@actionproducts.com

www.actionproducts.com



actionproducts.com/blogs



Product Warranty

Action® O.R. Patient Positioners are warranted to be free of defective materials and workmanship at the time of purchase and are serviceable for a period of two years provided the product is not subject to misuse, negligence, abuse, alteration, or failure to follow the Care and Use instructions (Instructions for Use). Accessories and disposable items are warranted to be free of defective materials and workmanship at the time of initial purchase only.

General Terms

- New accounts subject to credit approval
- Action accepts Visa®, MasterCard® and American Express®
- Handling, drop-ship, or other fees may apply
- Purchase order acceptance is confirmed by a sales order confirmation email. Customers are responsible for reviewing the order confirmation to verify details are accurate and acceptable. Contact Customer Service if the order confirmation is not received within two (2) business days or if there are errors in the sales order
- FCA – Action Products, Inc. INCOTERMS® 2020
- Specifications, prices, terms, and product availability are subject to change without notice
- Contracted terms may differ from General Terms

Product Returns

Please visit the [Product Return Info](#) page on our website for complete information.

- Action Products must pre-authorize all product returns by providing a required RMA number
- Items must be unused and returned in original product packaging
- A 20% restocking fee will be applied
- Customer is responsible for return shipping charges